



10/040791

COFC
HR

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

In re patent of:)
Geoffrey N. Roth et al.)
Patent No. 6,787,332) Group Art Unit:
Issued: 09/07/2004)
Title: Test Media & Quantitative or) Examiner: Henry Randall
Qualitative Method for Identification et al.) Cecelia Newman

REQUEST FOR RECONSIDERATION OF REFUSAL OF
CORRECTION OF PATENT CERTIFICATE

M/S Decisions & Certificates of Correction Branch
Commissioner for Patents
P.O. Box 1450
Alexandria, VA 22313-1450

Certificate
NOV 02 2005
of Correction


Dear Sirs:

The owner of U.S. Patent No. 6,787,332 for TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF *E. COLI*, GENERAL COLIFORMS, SALMONELLA, AND AEROMONAS IN A TEST SAMPLE, requests reconsideration of the refusal of correction of the patent and certificate of correction dated October 11, 2005.

Applicant has previously submitted three correct requests for correction. Attached hereto is the latest request, dated September 23, 2005, which includes all relevant information relating to this request, and provides the basis to show that the requested correction is supported by the prosecution history.

If the Examiner does not agree with the correction, Applicant requests that he contact the undersigned by telephone at 574-239-1966.

Respectfully submitted,

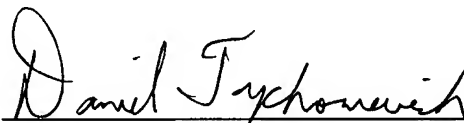

Daniel Tychonievich
Registration No. 41,358



DT:MLWH

Baker & Daniels LLP
Suite 250
205 West Jefferson Boulevard
South Bend, IN 46601
Telephone: 574-234-4149
Facsimile: 574-239-1900

Enclosure

CERTIFICATION OF MAILING

I hereby certify that this correspondence is being deposited with the United States Postal Service as First Class Mail in an envelope addressed to: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450, on: October 26, 2005.

DANIEL TYCHONIEVICH, REG. NO. 41,358

Name of Registered Representative

Daniel Tychonievich
Signature

October 26, 2005

Date



IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

In re patent of:)	
Geoffrey N. Roth et al.)	
Patent No. 6,787,332)	Group Art Unit:
Issued: 09/07/2004)	
Title: Test Media & Quantitative or)	Examiner: Henry Randall
Qualitative Method for Identification et al.)	Cecelia Newman

REQUEST FOR CORRECTION OF PATENT CERTIFICATE

M/S Decisions & Certificates of Correction Branch
Commissioner for Patents
P.O. Box 1450
Alexandria, VA 22313-1450

Dear Sirs:

The owner of U.S. Patent No. 6,787,332 for TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF *E. COLI*, GENERAL COLIFORMS, SALMONELLA, AND AEROMONAS IN A TEST SAMPLE, requests correction of the patent and certificate of correction dated May 24, 2005 as follows:

In Claim 3, Line 53, delete "3" and insert - β -.

In Claim 5, Line 60, delete " α " and insert - β -

In Claim 5, Line 63, delete " α " and insert - β -, such that claim 5 should now read:


5. The test medium of claim 72, wherein said first substrate consists essentially of 8-hydroxyquinoline- β -D-glucuronide, said second substrate consists essentially of 5-bromo-4-chloro-3-indole- α -D-galactoside, and said third substrate consists essentially of 6-chloro-3-indole- β -D-galactoside.

In connection therewith the Applicant offers the following information:

1. In a communication dated March 19, 2004, Applicant responded to an outstanding office action in which claim 80 of the application (now claim 5 of the patent), was amended to change its dependency (*see attached Exhibit A*). No other changes were made to claim 80.
2. On April 15, 2004, the Examiner issued a Notice of Allowance, allowing claims 72, 73, 76, 79-83 in response to the Amendment of March 19, 2004 (*see attached Exhibit B*). As can be seen the examiner did not make any changes to the text of any claims.
3. On September 7, 2004 the subject patent was issued under Patent No. 6,787,332 (*see attached Exhibit C*). As can be seen, Claim 5 in the patent contained an error on line 62 of the patent in which the Office had printed " β " where an α - should have been.
4. On March 2, 2005, Applicant requested a certificate of correction to rectify the above error (*see attached Exhibit D*). Applicant requested only that Line 62 of Claim 5 be corrected, and included the complete text of how claim 5 should read.
5. On May 24, 2005, Applicant received a Certificate of Correction for the subject patent, that indicated all occurrences of the " β " in claim 5 were changed to α - (*see attached Exhibit E*).
6. On June 6, 2005, Applicant requested a second Certificate of Correction to correct the Certificate issued on May 24, 2005 containing the erroneous replacement of all " β " to α -. Applicant also requested a correction of Claim 3, which, due to Applicant's error, contained a "3" in line 53, instead of the original β - (*see attached Exhibit F*).
7. On July 26, 2005 that request was denied stating that claims 3 and 5 were printed in accordance with the record (*see attached Exhibit G*). While Applicant will concede that claim 3 was printed in accordance with the amendment of March 19, 2004 (which contained the typographical error), it is clear from the record that claim 5 was printed in error in the actual patent and was later corrected in error as well.

Therefore, Applicant requests that a new certificate of correction be issued in accordance with the above request. Applicant hereby submits the fee of \$100 relating to the claim 3 correction. The office is hereby authorized to charge any additional fees or credit any overpayments to Deposit Account No. 02-0387, Baker & Daniels LLP. If the Examiner has any questions, he is invited to contact the undersigned at 574-234-4149.

Respectfully submitted,

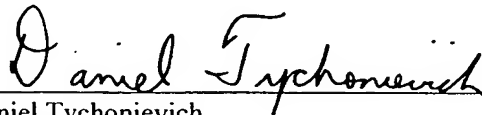

Daniel Tychonievich
Registration No. 41,358

DT:MLWH

Baker & Daniels LLP
Suite 250
205 West Jefferson Boulevard
South Bend, IN 46601
Telephone: 574-234-4149
Facsimile: 574-239-1900


Enc. Exhibits A - G
Check No. 104310 (\$100)

CERTIFICATION OF MAILING

I hereby certify that this correspondence is being deposited with the United States Postal Service as First Class Mail in an envelope addressed to:
Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450, on: September , 2005.

DANIEL TYCHONIEVICH, REG. NO. 41,358

Name of Registered Representative

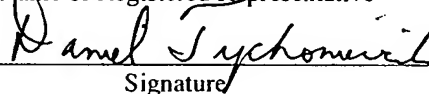

Signature

September 23, 2005

Date



The U. S. Patent and Trademark Office date-stamp hereon is acknowledgment that on the date stamped, the U. S. Patent and Trademark Office received:

TYPE OF PAPER	: Request for Correction of Patent Certificate
IN RE PATENT OF	: Geoffrey N. Roth <i>et al.</i>
PATENT NO.	: 6,787,332
ISSUED	: 09/07/2004
TITLE	: TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF <i>E. COLI</i> , GENERAL COLIFORMS, SALMONELLA, AND AEROMONAS IN A TEST SAMPLE
ENCLOSURES	: Request for Correction of Patent Certificate; Check No. 104310 in the amount of \$100.00; Exhibits A - G
DATE MAILED	: 09/23/2005
OUR REF.	: MLB0001.01 (65514.1)
DT/MLWH/mr	
BDDB01 3227327v1	




The U. S. Patent and Trademark Office date-stamp hereon is acknowledgment that on the date stamped, the U. S. Patent and Trademark Office received:

TYPE OF PAPER : Request for Correction of Patent Certificate
IN RE PATENT OF : Geoffrey N. Roth *et al.*
PATENT NO. : 6,787,332
ISSUED : 09/07/2004
TITLE : TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD
FOR IDENTIFICATION AND DIFFERENTIATION
OF E. COLI, GENERAL COLIFORMS, SALMONELLA, AND
AEROMONAS IN A TEST SAMPLE
ENCLOSURES : Request for Correction of Patent Certificate;
Check No. 104310 in the amount of \$100.00;
Exhibits A - G

DATE MAILED : 09/23/2005
OUR REF. : MLB0001.01 (65514.1)
DT/MLWH/mr

BDD01 3227327v1

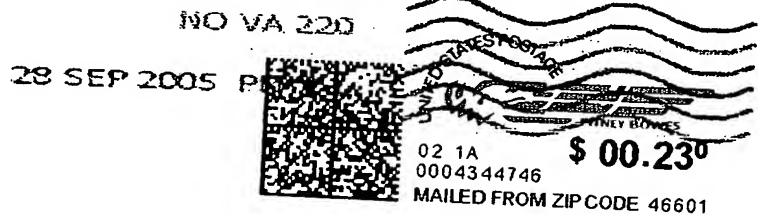


Baker & Daniels LLP
205 W. Jefferson, STE. 250
South Bend, IN 46601

RECEIVED
OCT 03 2005
BAKER & DANIELS

48801-1828



Application No. 10/040,791
Amended March 16, 2004
Reply to Office Action of November 28, 2003

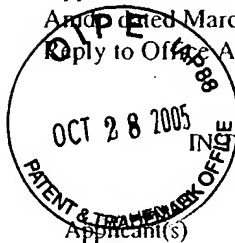

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Applicant(s) : Geoffrey N. Roth
Application No. : 10/040,791
Filed : January 7, 2002
Title : TEST MEDIA AND QUANTITATIVE OR QUALITATIVE
METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF
BIOLOGICAL MATERIALS IN A TEST SAMPLE
Group/Art Unit : 1651
Examiner : Ralph J. Gitomer
Docket No. : MLB0001.01

M/S FEE AMENDMENT
Commissioner for Patents
P.O. Box 1450
Alexandria, VA 22313-1450

AMENDMENT

Sir:

In response to the Office Action of October 22, 2003, please amend the above-identified application as follows.

Amendments to the Specification begins on page 2 of this paper.

Amendments to the Claims are reflected in the Listing of Claims which begins on page 3 of this paper.

Remarks/Arguments begin on page 5 of this paper.

AMENDMENTS TO THE SPECIFICATION

Please replace the title and add the following new paragraph before the paragraph beginning on line 1 of page 1:

Title: TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF E. COLI, GENERAL COLIFORMS, SALMONELLA, AND AEROMONAS IN A TEST SAMPLE

This Application claims priority from and is a continuation of U.S. Patent Application No. 09/357,606, which is incorporated in its entirety herein by reference, as filed on July 20, 1999 and subsequently issued as Patent No. 6,350,588.

LISTING OF CLAIMS

72. (Currently Amended) A test medium for detecting, quantifying and or differentiating general coliforms, *E. coli*, and ~~at least one of the genera *Aeromonas* and *Salmonella*~~, said test medium comprising: a nutrient base medium including ions of a salt; a first substrate which forms a first component of a first color in the presence of *E. coli*; a second substrate which forms a second component of a second color in the presence of ~~general coliforms~~ *Salmonella*; and a third substrate which forms a third component of a third color in the presence of ~~one of the genera *Aeromonas* and *Salmonella*~~; said second and third substrate forming said second and third components, respectively to make a fourth color, which is a combination of said second and third colors, in the presence of general coliforms, all of said colors being distinguishable from one another; said first substrate being a β -glucuronide nonchromogenic substrate; said second substrate being an α -D-galactoside chromogenic substrate and said third substrate[s] being a β -D-galactoside chromogenic substrate, ~~ehromogenic substrates~~; and said first color being substantially black.

73. (Original) The test medium of claim 72, wherein said first substrate is selected from the group consisting of 8-hydroxyquinoline- β -D-glucuronide, an esculetin glucuronide, and cyclohexenoesculetin- β -D-glucuronide.

74. (Cancelled)

75. (Cancelled)

76. (Currently Amended) The test medium of claim 72, wherein said first substrate is 8-hydroxyquinoline- β -D-glucuronide and forms a substantially nondiffusible compound in the presence of ions of said salt and *E. coli*, ~~and said third substrate also forms said third component of said third color in the presence *Shigella*~~.

77. (Cancelled)

78. (Cancelled)

79. (Original) The test medium of claim 72, wherein said salt comprises a metallic salt and said first component is water insoluble as formed by reaction with said ions.

80. (Currently Amended) The test medium of claim ~~75~~ 72, wherein said first substrate consists essentially of 8-hydroxyquinoline- β -D-glucuronide, said second substrate consists essentially of 5-bromo-4-chloro-3-indole- α -D-galactoside, and said third substrate consists essentially of 6-chloro-3-indole- β -D-galactoside.

81. (Currently Amended) A method for detecting, quantifying, and or differentiating colonies of *Aeromonas* from selected other biological entities in a test

sample, said method comprising the following steps: providing a base medium including ions of salt, a β -D-galactoside substrate that forms a first component of a first color in the presence of a first enzyme, an α -D-galactoside substrate that forms a second component of a second color distinguishable from said first color in the presence of a second enzyme, and a β -glucuronide nonchromogenic substrate that forms a third substantially black component in the presence of a third enzyme; inoculating the test medium with a test sample; incubating the test medium; and examining the test medium whereby aggregations of colonies of *Aeromonas* are indicated by said first color, and aggregations of colonies of *Salmonella* are indicated by said second color, and whereby colonies of general coliforms are indicated by a third color, said third color being a combination of said first and second colors.

82. (Original) The method as set forth in claim 81, further comprising the step of examining the test medium for *E. coli* as indicated by the presence of substantially black aggregates.

83. (Currently Amended) The method as set forth in Claim 82, wherein said β -glucuronide substrate is 8-hydroxyquinoline, ~~and further comprising the step of examining the test medium for *Shigella*, which are also indicated by said second color.~~

REMARKS

The Office Action of November 20, 2003 is acknowledged. Claims 72-83 were rejected in the Office Action.

Interview Summary

Applicant thanks the Examiner for participating in a telephone interview of February 12, 2004 in which the inventors Jonathan Roth and Geoffrey Roth were in attendance. Applicant provided photographs of the inventive test medium shown indicating *E.coli*, general coliforms, *Salmonella*, and *Aeromonas*, as four separate colors that are easily detected, differentiated and quantified. The Examiner concurred that Applicant's test medium was novel and non-obvious. However, the Examiner noted concerns with the conjunctives used in some of the claims and suggested modifications, which have been incorporated in the above amendments. Furthermore, the Examiner noted that the type of substrates should be identified in claim 72 to clarify the claim and clearly differentiate from the cited prior art. Accordingly, claim 72 has been amended to identify the type of substrates in the test medium.

Claim of Priority

The Examiner requested that the specification identify Application No. 09/357,606 filed July 7, 1999 from which subject Application claims priority. An appropriate amendment has been included.

Rejections Under 35 U.S.C. § 103(a)

Claims 72-75 and 77-83 were rejected under 35 U.S.C. § 103(a) as being obvious over the combination of U.S. Patent No. 5,726,231 to Roth *et al.* in view of the article to James *et al.* entitled "Detection of Specific Bacterial Enzymes by High Contrast Metal Chelate Formation: Specific Detection of *E. coli* on Multipoint Inoculated Plates Using 8-Hydroxyquinoline-Beta-D-Glucuronide" and in further view of U.S. Patent No. 4,308,348 to Monger, and PCT International Patent Application WO 98/55644 to Perry *et al.* As discussed in the interview, Roth *et al.* only teaches the use of two chromogens which act on *E. coli* and coliforms. A ph-indicator not a chromogen is used to identify the *Salmonella* and *Proteus*. Roth *et al.* does not teach or suggest the use of more than two enzyme specific substrates; otherwise, he would not have had to use a ph-indicator if he thought additional substrates would work. Furthermore, Roth *et al.* does not mention the detection of *Aeromonas*.

As for James *et al.*, James *et al.* only uses Hydroxyquinoline compound to attempt to identify *E. coli*. Furthermore, James *et al.* does not mention that Hydroxyquinoline can or should

be combined with chromogens or Indoxyl compounds. Furthermore, *Aeromonas* is not mentioned or suggested in any way in the *James*' article.

In addition, the test medium and method found in *Monget* is nothing like Applicant's test medium that includes two chromogenic and one non-chromogenic substrate in a nutrient based medium for growing bacteria and wherein bacteria can easily be detected, quantified and differentiated. In *Monget*, bacteria are grown separately in a gelose medium and then placed in a microtube (column 1, lines 50-52). *Monget* discloses a test medium where a β -galactosidase and β -glucosidase are impregnated on a disk of paper contained in the microtube. Although *Monget* mentions that the test may be also conducted with other substrates, *Monget* offers no teaching on how the bacteria would be differentiated in a test medium with more than two substrates. Furthermore, *Monget* offers no teaching on the use of a non-chromogenic substrate in combination with two chromogenic substrates in the test medium or a test medium for testing for *Aeromonas*.

Perry et al. discloses a test medium using α -D-galactosidase and β -D-galactosidase substrates for identifying the presence of *Salmonella*. The β -D-galactosidase substrate is cyclohexenoescluletin, which is a non-chromogenic substrate, although *Perry et al.* identifies it as a chromogenic substrate. *Perry's* medium cannot distinguish and differentiate *E.coli*, general coliforms and *Aeromonas* as all are positive to β -D-galactosidase and all would appear as substantially black colonies. Furthermore, *Perry et al.* offers no suggestion on how the substrates could be combined with another chromogenic substrate in order to differentiate or identify additional bacteria. In addition, *Perry et al.* does not teach the use of a β -D-glucuronide substrate nor a substrate for detecting *Aeromonas*. Accordingly, using the β -D-galactoside suggested by *Perry et al.* would render the test medium in Applicant's invention nonfunctional. As such, the above remarks and distinctions in clarifying claim amendments to claims 72, 73, 76 and 79-83 should be allowable over *Roth et al.* in view of *James et al.* and in further view of *Monget* and *Perry et al.*

Claim 76 was rejected under 35 U.S.C. § 103(a) as being obvious over a combination of *Roth* in view of *James* and further view of *Monget* and *Perry* as applied to claims 72-75, and 77-83 and in further view of *Kampfer*. Applicant notes that the rejection is now moot as the reference to detecting *Shigella* in claim 76 has been eliminated in the above amendment.

Claims 81-83 were rejected as being unpatentable over PCT International Publication No. WO 96/40980 to *Townsend et al.* in view of *Monget*. Applicant acknowledges that *Townsend et al.* mentions the detection of numerous bacteria including *Aeromonas* and also lists a multitude of various substrates in Table I. As discussed in the interview; however, *Townsend et al.* does not teach a method for detecting, quantifying or differentiating colonies of *Aeromonas* from other biological entities in a test sample. Although *Townsend et al.* mentions *Aeromonas*, he offers no explanation or teaching how *Aeromonas* should be quantified or differentiated from other bacteria.

Furthermore, *Townsend et al.* fails to teach the test medium used in the claimed method including a β -D-galactoside, α -D-galactoside and a β -D-glucuronide. Furthermore, in particular, claim 1 includes the method for detecting or differentiating *Aeromonas*, *Salmonella* and general coliforms. *Townsend et al.*, however, is directed towards merely quantifying the total number of viable bacteria present in ground beef and chicken (page 14, lines 21-23 and page 15, 19-22).

As discussed in the interview and above, *Monget* adds nothing to cure the lack of teachings in *Townsend et al.* regarding the subject method and test medium.

Claim Rejections Under 35 U.S.C. § 112, First Paragraph

Claims 72-83 were rejected as being non-enabling for other non-chromogenic substrates besides hydroxyquinoline- β -D-glucuronide. Applicant respectfully disagrees with this assertion as the use of other non-chromogenic substrates are specified in the specification on page 28, lines 11-17 and in Table V.

The Examiner also asserted that the three substrates required in the medium of the claims were not defined in the claims in a manner, which would permit one skilled in the art to make and use the invention. The substrates of the medium have been further defined and clarified in the attached amendments.

The Examiner also asserted that the claims are not enabled for distinguishing *Shigella* from *E. coli*. Applicant notes that this rejection is moot as the references to *Shigella* in claims 76 and 83 have been eliminated.

The Examiner also asserted that the specification is enabling for a method to detect and distinguish *E. coli*, coliforms and *Salmonella*, but does not reasonably provide an enablement for similar medium for detecting and distinguishing *E. coli*, coliforms and *Aeromonas*. As discussed in the interview and demonstrated in the photographs provided to the Examiner, Applicant's medium clearly provides a way of detecting and distinguishing *E. coli*, coliforms, *Aeromonas* and *Salmonella*. The photographs provided to the Examiner demonstrate the use of one example of the test medium as described on page 20, lines 20-25 in the Application. With this medium, *E. coli* shows as a substantially black color, general coliforms as a blue-violet color, *Aeromonas* as a pink color and *Salmonella* as a teal color as described on page 22, lines 14-31 and page 23, lines 1-8 of the specification. The Examiner is correct that *Aeromonas* does not produce β -glucuronide; however, as indicated in the specification on page 22, lines 26-31, *Aeromonas* is β -galactosidase positive, providing the pink color associated with 6-chloro-3-indolyl- β -D-galactoside.

Claim Rejections Under 35 U.S.C. § 112, Second Paragraph

Claim 74 was rejected for the redundant use of "one of". Applicant notes that this rejection is moot based on the cancellation of claim 74.

Claim 76 was rejected based on the assertion that "substantially nondefusable" is not defined by the claim or specification. Applicant respectfully traverses this rejection as the Examiner is aware the term "substantially" has been found to be definite as one of ordinary skill in the art knows what terms such as "substantially equal" mean. MPEP § 2173.05(b)(D).

The Examiner also objected that the title of the invention is not aptly descriptive. Accordingly, Applicant has provided a new title in the attached amendment.


Having fully responded to the Office Action, Applicant respectfully requests that remaining claims 72, 73, 76, and 79-83 be allowed. Applicant notes that certain claims have been amended and/or cancelled solely to advance prosecution of this Application and to obtain allowance on clearly allowable claims at the earliest possible date. Therefore, no admission may be inferred by the cancellations and amendments to the claims herein. Should the Examiner have any further questions or comments, he is invited to call the undersigned representative of the Applicant.



Respectfully submitted,



Daniel Tychonievich
Registration No. 41,358
Attorney for Applicant
BAKER & DANIELS
205 West Jefferson Blvd., Suite 250
South Bend, IN 46601
Tel: (574)234-4149
Fax: (574)239-1900

Enclosure: Postcard

I hereby certify that this correspondence is being deposited with the U.S. Postal Service as First Class Mail in an envelope addressed to: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450, on March <u>16</u> , 2004	
Date	<u>3/16/04</u>
 Daniel Tychonievich, Reg. No. 41,358	

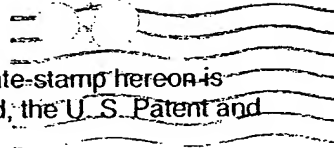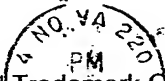



The U. S. Patent and Trademark Office date-stamp hereon is acknowledgment that on the date stamped, the U. S. Patent and Trademark Office received.

TYPE OF PAPER: Amendment
APPLICANT(S): Geoffrey N. Roth
TITLE: TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF BIOLOGICAL MATERIALS IN A TEST SAMPLE
SERIAL NO.: 10/040,791
FILING DATE: March 16, 2004
ENCLOSURES:

Check No. 102686 in the amount of \$55. (one month extension \$55.00)
Transmittal Form; Fee Transmittal for FY2004; Amendment (8 pages); Petition for Extension of Time Under 37 CFR 1.136(a)

Date Mailed: March 16, 2004
DT/kv
MLB0001.01 (65514.1)

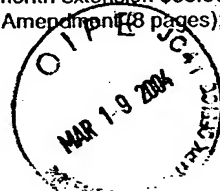

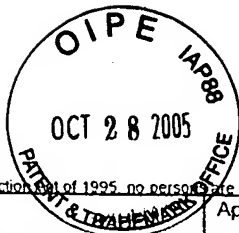

PTO/SB/21 (01-03)

Approved for use through 04/30/2003. OMB 0651-0031

U.S. Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE

Under the Paperwork Reduction Act of 1995, no person shall be required to respond to a collection of information unless it displays a valid OMB control number.

**TRANSMITTAL
FORM**

(to be used for all correspondence after initial filing)

Application Number	10/040,791
Filing Date	January 7, 2002
First Named Inventor	Geoffrey N. Roth
Art Unit	1651
Examiner Name	Ralph J. Gitomer
Attorney Docket Number	MLB0001.01

Total Number of Pages in This Submission

4

ENCLOSURES (Check all that apply)

<input checked="" type="checkbox"/> Fee Transmittal Form	<input type="checkbox"/> Drawing(s)	<input type="checkbox"/> After Allowance Communication to Group
<input checked="" type="checkbox"/> Fee Attached	<input type="checkbox"/> Licensing-related Papers	<input type="checkbox"/> Appeal Communication to Board of Appeals and Interferences
<input checked="" type="checkbox"/> Amendment/Reply	<input type="checkbox"/> Petition	<input type="checkbox"/> Appeal Communication to Group (Appeal Notice, Brief, Reply Brief)
<input type="checkbox"/> After Final	<input type="checkbox"/> Petition to Convert to a Provisional Application	<input type="checkbox"/> Proprietary Information
<input type="checkbox"/> Affidavits/declaration(s)	<input type="checkbox"/> Power of Attorney, Revocation Change of Correspondence Address	<input type="checkbox"/> Status Letter
<input checked="" type="checkbox"/> Extension of Time Request	<input type="checkbox"/> Terminal Disclaimer	<input checked="" type="checkbox"/> Other Enclosure(s) (please identify below):
<input type="checkbox"/> Express Abandonment Request	<input type="checkbox"/> Request for Refund	Return Postcard
<input type="checkbox"/> Information Disclosure Statement	<input type="checkbox"/> CD, Number of CD(s) _____	
<input type="checkbox"/> Certified Copy of Priority Document(s)	Remarks	
<input type="checkbox"/> Response to Missing Parts/Incomplete Application		
<input type="checkbox"/> Response to Missing Parts under 37 CFR 1.52 or 1.53		

SIGNATURE OF APPLICANT, ATTORNEY, OR AGENT


Firm or Individual	Daniel Tychonievich, Reg. No. 41,358 Baker & Daniels
Signature	<i>Daniel Tychonievich</i>
Date	March 16, 2004

CERTIFICATE OF TRANSMISSION/MAILING

I hereby certify that this correspondence is being facsimile transmitted to the USPTO or deposited with the United States Postal Service with sufficient postage as first class mail in an envelope addressed to: Commissioner for Patents, Washington, DC 20231 on this date: March 16, 2004			
Typed or printed	Daniel Tychonievich		
Signature	<i>Daniel Tychonievich</i>	Date	March 16, 2004

This collection of information is required by 37 CFR 1.5. The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 12 minutes to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, Washington, DC 20231. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Commissioner for Patents, Washington, DC 20231.

If you need assistance in completing the form, call 1-800-PTO-9199 (1-800-786-9199) and select option 2.

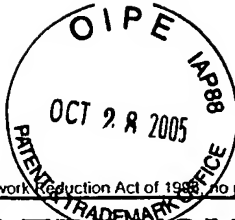

PTO/SB/17 (10-03)

Approved for use through 07/31/2006. OMB 0651-0032

U.S. Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE

Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it displays a valid OMB control number.

FEE TRANSMITTAL for FY 2004

Effective 10/01/2003. Patent fees are subject to annual revision.

☒ Applicant claims small entity status. See 37 CFR 1.27

TOTAL AMOUNT OF PAYMENT (\$ 55.00

Complete if Known

Application Number	10/040,791
Filing Date	January 7, 2002
First Named Inventor	Geoffrey N. Roth
Examiner Name	Ralph J. Gitomer
Art Unit	1651
Attorney Docket No.	MLB0001.01

METHOD OF PAYMENT (check all that apply)☒ Check ☐ Credit card ☐ Money Order ☐ Other ☐ None☒ Deposit Account:Deposit Account Number
02-0387
Deposit Account Name
Baker & Daniels

The Director is authorized to: (check all that apply)

☒ Charge fee(s) indicated below ☒ Credit any overpayments☐ Charge any additional fee(s) or any underpayment of fee(s)☐ Charge fee(s) indicated below, except for the filing fee to the above-identified deposit account.**FEE CALCULATION****1. BASIC FILING FEE**

Large Entity Fee Code (\$)	Small Entity Fee Code (\$)	Fee Description	Fee Paid
1001 770	2001 385	Utility filing fee	
1002 340	2002 170	Design filing fee	
1003 530	2003 265	Plant filing fee	
1004 770	2004 385	Reissue filing fee	
1005 160	2005 80	Provisional filing fee	

SUBTOTAL (1) (\$)

2. EXTRA CLAIM FEES FOR UTILITY AND REISSUE

Total Claims	Extra Claims	Fee from below	Fee Paid
Independent Claims	- 20** =	X	
Multiple Dependent Claims	- 3** =	X	

Large Entity Fee Code (\$)	Small Entity Fee Code (\$)	Fee Description
1202 18	2202 9	Claims in excess of 20
1201 86	2201 43	Independent claims in excess of 3
1203 290	2203 145	Multiple dependent claim, if not paid
1204 86	2204 43	** Reissue independent claims over original patent
1205 18	2205 9	** Reissue claims in excess of 20 and over original patent

SUBTOTAL (2) (\$)

**or number previously paid, if greater. For Reissues, see above

FEE CALCULATION (continued)**3. ADDITIONAL FEES**

Large Entity Small Entity

Fee Code (\$)	Fee Code (\$)	Fee Description	Fee Paid
1051 130	2051 65	Surcharge - late filing fee or oath	
1052 50	2052 25	Surcharge - late provisional filing fee or cover sheet	
1053 130	1053 130	Non-English specification	
1812 2,520	1812 2,520	For filing a request for ex parte reexamination	
1804 920*	1804 920*	Requesting publication of SIR prior to Examiner action	
1805 1,840*	1805 1,840*	Requesting publication of SIR after Examiner action	
1251 110	2251 55	Extension for reply within first month	55.00
1252 420	2252 210	Extension for reply within second month	
1253 950	2253 475	Extension for reply within third month	
1254 1,480	2254 740	Extension for reply within fourth month	
1255 2,010	2255 1,005	Extension for reply within fifth month	
1401 330	2401 165	Notice of Appeal	
1402 330	2402 165	Filing a brief in support of an appeal	
1403 290	2403 145	Request for oral hearing	
1451 1,510	1451 1,510	Petition to institute a public use proceeding	
1452 110	2452 55	Petition to revive - unavoidable	
1453 1,330	2453 665	Petition to revive - unintentional	
1501 1,330	2501 665	Utility issue fee (or reissue)	
1502 480	2502 240	Design issue fee	
1503 640	2503 320	Plant issue fee	
1460 130	1460 130	Petitions to the Commissioner	
1807 50	1807 50	Processing fee under 37 CFR 1.17(q)	
1806 180	1806 180	Submission of Information Disclosure Stmt	
8021 40	8021 40	Recording each patent assignment per property (times number of properties)	
1809 770	2809 385	Filing a submission after final rejection (37 CFR 1.129(a))	
1810 770	2810 385	For each additional invention to be examined (37 CFR 1.129(b))	
1801 770	2801 385	Request for Continued Examination (RCE)	
1802 900	1802 900	Request for expedited examination of a design application	

Other fee (specify)

*Reduced by Basic Filing Fee Paid

SUBTOTAL (3) (\$ 55.00


SUBMITTED BY

Name (Print/Type)	Daniel Tychonievich	Registration No. (Attorney/Agent)	41,358	Telephone	574/234-4149
Signature	<i>Daniel Tychonievich</i>	Date	March 16, 2004		

WARNING: Information on this form may become public. Credit card information should not be included on this form. Provide credit card information and authorization on PTO-2038.

This collection of information is required by 37 CFR 1.17 and 1.27. The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 12 minutes to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, P.O. Box 1450, Alexandria, VA 22313-1450. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2

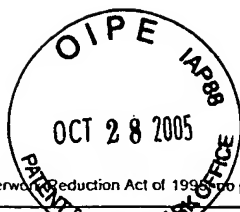

Under the paperwork Reduction Act of 1995 no persons are required to respond to a collection of information unless it displays a valid OMB control number.

PTO/SB/22 (08-03)

Approved for use through 7/31/2006. OMB 0651-0031

U.S. Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE

PETITION FOR EXTENSION OF TIME UNDER 37 CFR 1.136(a)		Docket Number (Optional) MLB0001.01	
In re Application of Geoffrey N. Roth			
Application Number 10/040,791		Filed January 7, 2002	
For TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF BIOLOGICAL MATERIALS IN A TEST SAMPLE			
Art Unit 1651		Examiner Ralph J. Gitomer	
<p>This is a request under the provisions of 37 CFR 1.136(a) to extend the period for filing a reply in the above identified application.</p> <p>The requested extension and appropriate non-small-entity fee are as follows (check time period desired):</p> <p><input checked="" type="checkbox"/> One month (37 CFR 1.17(a)(1)) \$ <u>110.00</u></p> <p><input type="checkbox"/> Two months (37 CFR 1.17(a)(2)) \$ _____</p> <p><input type="checkbox"/> Three months (37 CFR 1.17(a)(3)) \$ _____</p> <p><input type="checkbox"/> Four months (37 CFR 1.17(a)(4)) \$ _____</p> <p><input type="checkbox"/> Five months (37 CFR 1.17(a)(5)) \$ _____</p> <p><input checked="" type="checkbox"/> Applicant claims small entity status. See 37 CFR 1.27. Therefore, the fee amount shown above is reduced by one-half, and the resulting fee is: \$ <u>55.00</u></p> <p><input checked="" type="checkbox"/> A check in the amount of the fee is enclosed.</p> <p><input type="checkbox"/> Payment by credit card. Form PTO-2038 is attached.</p> <p><input type="checkbox"/> The Director has already been authorized to charge fees in this application to a Deposit Account.</p> <p><input checked="" type="checkbox"/> The Director is hereby authorized to charge any fees which may be required, or credit any overpayment, to Deposit Account Number <u>02-0387</u></p> <p>I have enclosed a duplicate copy of this sheet.</p> <p>I am the <input type="checkbox"/> applicant/inventor.</p> <p><input type="checkbox"/> assignee of record of the entire interest. See 37 CFR 3.71. Statement under 37 CFR 3.73(b) is enclosed (Form PTO/SB/96).</p> <p><input checked="" type="checkbox"/> attorney or agent of record. Registration Number <u>41,358</u></p> <p><input type="checkbox"/> attorney or agent under 37 CFR 1.34(a). Registration number if acting under 37 CFR 1.34(a) _____</p> <p>WARNING: Information on this form may become public. Credit card information should not be included on this form. Provide credit card information and authorization on PTO-2038.</p> <p><u>March 16, 2004</u> Date</p> <p><u>(574) 234-4149</u> Telephone Number</p> <p><u>Daniel Tychonievich</u> Signature</p> <p><u>Daniel Tychonievich</u> Typed or printed name</p> <p>NOTE: Signatures of all the inventors or assignees of record of the entire interest or their representative(s) are required. Submit multiple forms if more than one signature is required, see below.</p> <p><input type="checkbox"/> Total of _____ forms are submitted.</p>			

This collection of information is required by 37 CFR 1.136(a). The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 6 minutes to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, P.O. Box 1450, Alexandria, VA 22313-1450. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2.



UNITED STATES PATENT AND TRADEMARK OFFICE

UNITED STATES DEPARTMENT OF COMMERCE
United States Patent and Trademark Office
Address: COMMISSIONER FOR PATENTS
P.O. Box 1450
Alexandria, Virginia 22313-1450
www.uspto.gov

NOTICE OF ALLOWANCE AND FEE(S) DUE

27187 7590 04/15/2004
BAKER & DANIELS
205 W. JEFFERSON BOULEVARD
SUITE 250
SOUTH BEND, IN 46601

RECEIVED

APR 19 2004

BAKER & DANIELS

EXAMINER

GITOMER, RALPH J

ART UNIT

PAPER NUMBER

1651

DATE MAILED: 04/15/2004

APPLICATION NO.	FILING DATE	FIRST NAMED INVENTOR	ATTORNEY DOCKET NO.	CONFIRMATION NO.
10/040,791	01/07/2002	Geoffrey N. Roth	MLB0001.01	3954

TITLE OF INVENTION: TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF BIOLOGICAL MATERIALS IN A TEST SAMPLE

APPLN. TYPE	SMALL ENTITY	ISSUE FEE	PUBLICATION FEE	TOTAL FEE(S) DUE	DATE DUE
nonprovisional	YES	\$665	\$300	\$965	07/15/2004

THE APPLICATION IDENTIFIED ABOVE HAS BEEN EXAMINED AND IS ALLOWED FOR ISSUANCE AS A PATENT. PROSECUTION ON THE MERITS IS CLOSED. THIS NOTICE OF ALLOWANCE IS NOT A GRANT OF PATENT RIGHTS. THIS APPLICATION IS SUBJECT TO WITHDRAWAL FROM ISSUE AT THE INITIATIVE OF THE OFFICE OR UPON PETITION BY THE APPLICANT. SEE 37 CFR 1.313 AND MPEP 1308.

THE ISSUE FEE AND PUBLICATION FEE (IF REQUIRED) MUST BE PAID WITHIN THREE MONTHS FROM THE MAILING DATE OF THIS NOTICE OR THIS APPLICATION SHALL BE REGARDED AS ABANDONED. THIS STATUTORY PERIOD CANNOT BE EXTENDED. SEE 35 U.S.C. 151. THE ISSUE FEE DUE INDICATED ABOVE REFLECTS A CREDIT FOR ANY PREVIOUSLY PAID ISSUE FEE APPLIED IN THIS APPLICATION. THE PTOL-85B (OR AN EQUIVALENT) MUST BE RETURNED WITHIN THIS PERIOD EVEN IF NO FEE IS DUE OR THE APPLICATION WILL BE REGARDED AS ABANDONED.

HOW TO REPLY TO THIS NOTICE:

I. Review the SMALL ENTITY status shown above.

If the SMALL ENTITY is shown as YES, verify your current SMALL ENTITY status:

A. If the status is the same, pay the TOTAL FEE(S) DUE shown above.

B. If the status is changed, pay the PUBLICATION FEE (if required) and twice the amount of the ISSUE FEE shown above and notify the United States Patent and Trademark Office of the change in status, or

If the SMALL ENTITY is shown as NO:

A. Pay TOTAL FEE(S) DUE shown above, or

B. If applicant claimed SMALL ENTITY status before, or is now claiming SMALL ENTITY status, check the box below and enclose the PUBLICATION FEE and 1/2 the ISSUE FEE shown above.

☐ Applicant claims SMALL ENTITY status.
See 37 CFR 1.27.

II. PART B - FEE(S) TRANSMITTAL should be completed and returned to the United States Patent and Trademark Office (USPTO) with your ISSUE FEE and PUBLICATION FEE (if required). Even if the fee(s) have already been paid, Part B - Fee(s) Transmittal should be completed and returned. If you are charging the fee(s) to your deposit account, section "4b" of Part B - Fee(s) Transmittal should be completed and an extra copy of the form should be submitted.

III. All communications regarding this application must give the application number. Please direct all communications prior to issuance to Mail Stop ISSUE FEE unless advised to the contrary.

IMPORTANT REMINDER: Utility patents issuing on applications filed on or after Dec. 12, 1980 may require payment of maintenance fees. It is patentee's responsibility to ensure timely payment of maintenance fees when due.



Notice of Allowability

Application No.

10/040,791

Examiner

Ralph Gitomer

Applicant(s)

ROTH ET AL.

Art Unit

1651

-- The MAILING DATE of this communication appears on the cover sheet with the correspondence address--
aims being allowable, PROSECUTION ON THE MERITS IS (OR REMAINS) CLOSED in this application. If not included
with (or previously mailed), a Notice of Allowance (PTOL-85) or other appropriate communication will be mailed in due course. THIS
ICE OF ALLOWABILITY IS NOT A GRANT OF PATENT RIGHTS. This application is subject to withdrawal from issue at the initiative
e Office or upon petition by the applicant. See 37 CFR 1.313 and MPEP 1308.

☒ This communication is responsive to amendment of 3/19/04.

☒ The allowed claim(s) is/are 72 73 76 79-83.

☐ The drawings filed on _____ are accepted by the Examiner.

☐ Acknowledgment is made of a claim for foreign priority under 35 U.S.C. § 119(a)-(d) or (f).

a) ☐ All b) ☐ Some* c) ☐ None of the:

1. ☐ Certified copies of the priority documents have been received.

2. ☐ Certified copies of the priority documents have been received in Application No. _____.

3. ☐ Copies of the certified copies of the priority documents have been received in this national stage application from the
International Bureau (PCT Rule 17.2(a)).

* Certified copies not received: _____.

plicant has THREE MONTHS FROM THE "MAILING DATE" of this communication to file a reply complying with the requirements
ted below. Failure to timely comply will result in ABANDONMENT of this application.
IIS THREE-MONTH PERIOD IS NOT EXTENDABLE.

☐ A SUBSTITUTE OATH OR DECLARATION must be submitted. Note the attached EXAMINER'S AMENDMENT or NOTICE OF
INFORMAL PATENT APPLICATION (PTO-152) which gives reason(s) why the oath or declaration is deficient.

☐ CORRECTED DRAWINGS (as "replacement sheets") must be submitted.

(a) ☐ including changes required by the Notice of Draftperson's Patent Drawing Review (PTO-948) attached

1) ☐ hereto or 2) ☐ to Paper No./Mail Date _____.

(b) ☐ including changes required by the attached Examiner's Amendment / Comment or in the Office action of
Paper No./Mail Date _____.

identifying indicia such as the application number (see 37 CFR 1.84(c)) should be written on the drawings in the front (not the back) of
each sheet. Replacement sheet(s) should be labeled as such in the header according to 37 CFR 1.121(d).

☐ DEPOSIT OF and/or INFORMATION about the deposit of BIOLOGICAL MATERIAL must be submitted. Note the
attached Examiner's comment regarding REQUIREMENT FOR THE DEPOSIT OF BIOLOGICAL MATERIAL.

achment(s)

☐ Notice of References Cited (PTO-892)

☐ Notice of Draftperson's Patent Drawing Review (PTO-948)

☐ Information Disclosure Statements (PTO-1449 or PTO/SB/08),
Paper No./Mail Date _____

☐ Examiner's Comment Regarding Requirement for Deposit
of Biological Material

5. ☐ Notice of Informal Patent Application (PTO-152)

6. ☒ Interview Summary (PTO-413),
Paper No./Mail Date _____

7. ☐ Examiner's Amendment/Comment

8. ☐ Examiner's Statement of Reasons for Allowance

9. ☐ Other _____

Ralph Gitomer
Ralph Gitomer
Primary Examiner
Art Unit: 1651



UNITED STATES PATENT AND TRADEMARK OFFICE

UNITED STATES DEPARTMENT OF COMMERCE
United States Patent and Trademark Office
Address: COMMISSIONER FOR PATENTS
P.O. Box 1450
Alexandria, Virginia 22313-1450
www.uspto.gov

APPLICATION NO.	FILING DATE	FIRST NAMED INVENTOR	ATTORNEY DOCKET NO.	CONFIRMATION NO.
10/040,791	01/07/2002	Geoffrey N. Roth	MLB0001.01	3954
27187	7590	04/15/2004	EXAMINER GITOMER, RALPH J	
BAKER & DANIELS 205 W. JEFFERSON BOULEVARD SUITE 250 SOUTH BEND, IN 46601			ART UNIT 1651	PAPER NUMBER

DATE MAILED: 04/15/2004

Determination of Patent Term Adjustment under 35 U.S.C. 154 (b)
(application filed on or after May 29, 2000)

The Patent Term Adjustment to date is 246 day(s). If the issue fee is paid on the date that is three months after the mailing date of this notice and the patent issues on the Tuesday before the date that is 28 weeks (six and a half months) after the mailing date of this notice, the Patent Term Adjustment will be 246 day(s).

If a Continued Prosecution Application (CPA) was filed in the above-identified application, the filing date that determines Patent Term Adjustment is the filing date of the most recent CPA.

Applicant will be able to obtain more detailed information by accessing the Patent Application Information Retrieval (PAIR) system (<http://pair.uspto.gov>).

Any questions regarding the Patent Term Extension or Adjustment determination should be directed to the Office of Patent Legal Administration at (703) 305-1383. Questions relating to issue and publication fee payments should be directed to the Customer Service Center of the Office of Patent Publication at (703) 305-8283.

Issue Fee Due
REPLY DUE: 7.15.04
DOCKET 1 KV **DOCKET 2** ml



US006787332B2

(12) **United States Patent**
Roth et al.

(10) **Patent No.:** **US 6,787,332 B2**
(45) Date of Patent: **Sep. 7, 2004**

(54) **TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR IDENTIFICATION AND DIFFERENTIATION OF *E. COLI*, GENERAL COLIFORMS, SALMONELLA, AND AEROMONAS IN A TEST SAMPLE**

(75) **Inventors:** **Geoffrey N. Roth, Goshen, IN (US); Jonathan N. Roth, Goshen, IN (US)**

(73) **Assignee:** **Micrology Laboratories, LLC, Goshen, IN (US)**

(*) **Notice:** Subject to any disclaimer, the term of this patent is extended or adjusted under 35 U.S.C. 154(b) by 246 days.

(21) **Appl. No.:** **10/040,791**

(22) **Filed:** **Jan. 7, 2002**

(65) **Prior Publication Data**

US 2002/0090668 A1 Jul. 11, 2002

Related U.S. Application Data

(63) Continuation of application No. 09/357,606, filed on Jul. 20, 1999, now Pat. No. 6,350,588.

(51) **Int. Cl.⁷** **C12Q 1/04**

(52) **U.S. Cl.** **435/34; 435/38**

(58) **Field of Search** **435/34, 38, 7.32, 435/7.35, 252.4, 252.8, 14**

(56) **References Cited**

U.S. PATENT DOCUMENTS

4,308,348 A * 12/1981 Monget 435/38
 5,726,031 A * 3/1998 Roth et al. 435/34
 6,350,588 B1 * 2/2002 Roth et al. 435/34

FOREIGN PATENT DOCUMENTS

WO WO 96/40980 * 12/1996
 WO WO 98/55644 * 12/1998

OTHER PUBLICATIONS

Sartory D. A Simple Membrane Enumeration Medium for Coliforms and *E. coli* Utilizing a Chromogenic Glucuronide Substrate. Proceedings 1992 Water Quality Technology Conference, 1992.*

McDaniels, A. Confirmational Identification of *E. coli*, ... Applied and Environmental Micro 1996 62(9)3350-3354.*

Kampfer P. Glycosidase Profiles of Members of the Family Enterobacteriaceae. J of Clinical Micro 1991 29(12)2877-2879.*

James A. Detection of Specific Bacterial Enzymes by High Contrast Metal Chelate Formation. Zbl. Bakr. Hyg A 267, 316-321, 1988.*

* cited by examiner

Primary Examiner—Ralph Gitomer

(74) *Attorney, Agent, or Firm*—Baker & Daniels

(57) **ABSTRACT**

A test medium and method for detecting, quantifying, identifying and differentiating up to four (4) separate biological materials in a test sample. A test medium is disclosed which allows quantifying and differentiating under ambient light aggregates of biological entities producing specific enzymes, which might include general coliforms, *E. coli*, *Aeromonas*, and *Salmonella* or *Shigella* in a single test medium. A new class of nonchromogenic substrate is disclosed which produce a substantially black, non-diffusible precipitate. This precipitate is not subject to interference from other chromogenic substrates present in the test medium. In a preferred form, the substrates are selected such that *E. coli* colonies present in the test medium show as substantially black, general coliforms colonies show in the test medium as a blue-violet color, *Aeromonas* colonies present in the test medium show as a generally red-pink color, and *Salmonella* or *Shigella* colonies show as a generally teal-green color. Other microorganisms and color possibilities for detection and quantification thereof are also disclosed. An inhibitor and method for making a test medium incorporating the inhibitor are disclosed.

8 Claims, No Drawings

1

**TEST MEDIA AND QUANTITATIVE OR
QUALITATIVE METHOD FOR
IDENTIFICATION AND DIFFERENTIATION
OF *E. COLI*, GENERAL COLIFORMS,
SALMONELLA, AND AEROMONAS IN A
TEST SAMPLE**

This Application claims priority from and is a continuation of U.S. patent application Ser. No. 09/357,606, which is incorporated in its entirety herein by reference, as filed on Jul. 20, 1999 and subsequently issued as U.S. Pat. No. 6,350,588.

BACKGROUND OF THE INVENTION

The present invention relates to a test medium and method for the detection, quantification, identification and/or differentiation of biological materials in a sample which may contain a plurality of different biological materials.

Bacteria are the causative factor in many diseases of humans, higher animals and plants, and are commonly transmitted by carriers such as water, beverages, food and other organisms. The testing of these potential carriers of bacteria is of critical importance and generally relies on "indicator organisms." Borrego et al., *Microbiol. Sem.* 13:413-426, (1998). For example, *Escherichia coli* (*E. coli*) is a gram negative member of the family Enterobacteriaceae which is part of the normal intestinal flora of warm blooded animals, and its presence indicates fecal contamination (e.g., raw sewage). Even though most strains of *E. coli* are not the actual cause of disease, their presence is a strong indication of the possible presence of pathogens associated with intestinal disease, such as cholera, dysentery, and hepatitis, among others. Consequently, *E. coli* has become a prime indicator organism for fecal contamination, and as a result, any method which differentiates and identifies *E. coli* from other bacteria is very useful.

Others members of the family Enterobacteriaceae, commonly referred to as "general coliforms," especially the genera *Citrobacter*, *Enterobacter* and *Klebsiella*, are also considered to be significant indicator organisms for the quality of water, beverages and foods. Therefore, tests to identify and differentiate general coliforms from *E. coli* are also very useful. Also, various species of the genus *Aeromonas* have been shown to not only be potential pathogens, but to have a correlation to other indicator organisms (Pettibone et al., *J Appl. Microbiol.* 85:723-730 (1998)). Current test methods to identify, separate and enumerate *Aeromonas* spp. from the very similar Enterobacteriaceae have been lacking and most of the current methods utilizing enzyme substrates do not separate *Aeromonas* spp. from Enterobacteriaceae due to their almost identical biochemical profiles. Any method that depends upon the identification of general coliforms by means of a β -galactosidase substrate either does not differentiate *Aeromonas* spp. from general coliforms or eliminates *Aeromonas* from the sample by the use of specific inhibitors (antibiotic such as cefsulodin). Brenner et al., *Appl. Envir. Microbio.* 59:3534-44 (1993). They do not differentiate, identify and enumerate *Aeromonas* along with *E. coli* and general coliforms. Landre et al., *Letters Appl. Microbiol.* 26:352-354 (1998). Improved test methods to effectively identify, separate and enumerate such

2

bacterial types are needed, and there is a continuing search for faster, more accurate, easier to use and more versatile test methods and apparatus in this area.

Numerous test methods have been utilized to determine, identify and enumerate one or more indicator organisms. Some of these test methods only indicate the presence or absence of the microorganism, while others also attempt to quantify one or more of the particular organisms in the test sample. For example, a qualitative test referred to as the Presence/Absence (or P/A) test, may be utilized to determine the presence or absence of coliforms and *E. coli* in a test sample. A test medium including the β -galactosidase substrate O-nitrophenyl- β -D-galactopyranoside (ONPG), and the β -glucuronidase substrate 4-methyl-umbrelliferyl- β -D-glucuronide (MUG), is inoculated with the test sample. To differentiate the general coliforms from *E. coli*, this test relies on the fact that generally all coliforms produce P-galactosidase, whereas only *E. coli* also produces β -glucuronidase in addition to β -galactosidase. If any coliforms are present (including *E. coli*), the broth medium turns a yellow color due to the activity of the galactosidase enzyme on the ONPG material, causing the release of a diffusible yellow pigment. If *E. coli* is present, the broth medium will demonstrate a blue fluorescence when irradiated with ultraviolet rays, due to the breakdown of the MUG reagent with the release of the fluorogenic dye caused by the production of the glucuronidase enzyme. These reactions are very specific, and allow the presence of both coliforms in general, as well as *E. coli* to be identified in a single sample. A disadvantage of this test is that it is not directly quantitative for either bacterial type, since both reagents produce diffusible pigments. A second disadvantage is that there may be a false positive coliform reaction if *Aeromonas* spp. are present in the test sample. This has been shown to be possible even when there are inhibitors present to supposedly prevent this from occurring (Landre et al., *Letters Appl. Microbiol.* 26:352-354 (1998)). The test also requires specific equipment for producing the ultraviolet rays. Further, this test may only be used to detect coliforms and *E. coli*. Other important microorganisms, such as the strain *E. coli* 0157 which is glucuronidase negative, are not detected, nor are other non-galactosidase-glucuronidase producing microorganisms.

The Violet Red Bile Agar (VRBA) method has been used to determine the quantity of both coliform and *E. coli* in a test sample. The test medium used in this method includes bile salts (to inhibit non-coliforms), lactose and the pH indicator neutral red. As coliforms (including *E. coli*) grow in the medium, the lactose is fermented with acid production, and the neutral red in the area of the bacterial colony becomes a brick red color. The results of this test are not always easy to interpret, and in order to determine the presence of *E. coli*, confirming follow-up tests, such as brilliant green lactose broth fermentation, growth in EC broth at 44.5° C. and streaking on Eosin Methylene Blue Agar (EMBBA), must be performed.

The Membrane Filter (MF) method utilizes micropore filters through which samples are passed so that the bacteria are retained on the surface of the filter. This method is used most often when bacterial populations are very small, and a large sample is needed to get adequate numbers. The filter is

then placed on the surface of a chosen medium, incubated, and the bacterial colonies growing on the membrane filter surface are counted and evaluated. This method is widely used and provides good results when combined with proper reagents and media. A disadvantage of this method is that it is expensive and time-consuming. It also does not work well with solid samples, or samples with high particulate counts. The MF method can be used in conjunction with the inventive method described in this application.

The m-Endo method is also used to determine the quantity of *E. coli* and general coliforms and is an official USEPA approved method for testing water quality. The medium is commonly used with a membrane filter and *E. coli* and general coliform colony forming units (CFU) grow as dark colonies with a golden green metallic sheen. Due to a proven high rate of false positive error, typical colonies must be confirmed by additional testing. *Standard Methods for the Examination of water and Wastewater*, 20th Edition, 9-10 & 9-60 (1998).

Other tests, such as the Most Probable Number (MPN), utilize lactose containing broths (LST, BGLB, EC) to estimate numbers of general coliforms and *E. coli*, but have also been shown to have high rates of error as well as being cumbersome and slow to produce results. Evans et al., *Appl. Envir. Microbiol.* 41:130-138 (1981).

The reagent 5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside (X-gal) is a known test compound for identifying coliforms. When acted on by the β -galactosidase enzyme produced by coliforms, X-gal forms an insoluble indigo blue precipitate. X-gal can be incorporated into a nutrient medium such as an agar plate, and if a sample containing coliforms is present, the coliforms will grow as indigo blue colonies. X-gal has the advantage over the compound ONPG, described above, in that it forms a water insoluble precipitate rather than a diffusible compound, thereby enabling a quantitative determination of coliforms to be made when the test sample is incorporated into or onto a solidified medium, or when coliform colonies grow on the surface of a membrane filter resting on a pad saturated with a liquid medium or on a membrane filter resting on a solid medium. Further, it does not require the use of ultraviolet light.

A similar compound, 5-bromo-4-chloro-3-indolyl- β -D-glucuronide (X-gluc) is a known test compound for identifying *E. coli*. When acted on by the β -glucuronidase enzyme produced by most *E. coli*, X-gluc forms an insoluble indigo blue precipitate. X-gluc has the advantage over the compound MUG, described above, in that it forms a water insoluble precipitate, rather than a diffusible compound, thereby enabling a quantitative determination of *E. coli* to be made when the test sample is incorporated into or onto a solidified medium. X-gluc and its ability to identify *E. coli* are described in Watkins, et al., *Appl. Envir. Microbiol.* 54:1874-1875 (1988). A similar compound, indoxyl- β -D-glucuronide, which also produces sharp blue colonies of *E. coli*, was described in Ley, et al., *Can. J. Microbiol.* 34:690-693 (1987).

Although X-gal and X-gluc are each separately useful in the quantitative determination of either coliforms (X-gal) or *E. coli* (X-gluc), these indicator compounds have the disadvantage that they each contain the same chromogenic

component. Therefore, they cannot be used together to identify and distinguish both *E. coli* and general coliforms in a single test with a single sample, since they both generate identically hued indigo blue colonies. A person using both reagents together would be able to quantitatively identify the total number of coliforms, the same as if X-gal were used alone, but would not be able to indicate which of the colonies were *E. coli* and which were other coliforms besides *E. coli*.

A recently developed and highly commercially successful test method and test medium for quantitatively identifying and differentiating general coliforms and *E. coli* in a test sample is described in U.S. Pat. Nos. 5,210,022, and 5,393,662, both of which share common inventorship with the present application and which are hereby incorporated by reference. This method and test medium improves upon prior art methods by allowing the quantitative identification of general coliforms and *E. coli* in a single sample. Additional confirmatory tests are not necessary. The test sample is added to a medium containing a β -galactosidase substrate, such as 6-chloroindolyl- β -D-galactoside, and a β -glucuronidase substrate, such as 5-bromo-4-chloro-3-indolyl- β -D-glucuronide (X-gluc). The β -galactosidase substrate is capable of forming a water insoluble precipitate of a first color upon reacting with β -galactosidase, and the β -glucuronidase substrate is capable of forming a water insoluble precipitate of a second color, contrasting with the first color, upon reacting with β -glucuronidase. As a result, general coliforms may be quantified by enumerating the colonies of the first color (having β -galactosidase activity), and *E. coli* may be quantified by enumerating the colonies of the second color (having both β -galactosidase and β -glucuronidase activity). This technology has been widely copied.

Another recently developed test method and apparatus provides excellent results for the differentiation and identification of general coliforms, *E. coli* and *E. coli* 0157 strains and non-coliform Enterobacteriaceae. The method and test medium are described in U.S. Pat. No. 5,726,031, which shares common inventorship with the present application, and which is hereby incorporated by reference.

A certain class of substrates, referred to herein as "nonchromogenic," have been used to detect various microorganisms. A dipslide for detecting *E. coli* using hydroxyquinoline- β -D-glucuronide is disclosed by Dalet et al., *J. Clin. Microbiol.* 33(5):1395-8 (1995). Similarly, a technique for detection of *E. coli* in an agar-based medium using 8-hydroxyquinoline- β -D-glucuronide is disclosed by James et al., *Zentralbl. Bakteriol. Mikrobiol. Hyg. [A]*, 267(3):316-21 (1988).

It is desirable to further improve the distinguishing colors generated in the test media. That is to say, in prior art test media which detect and distinguishing two biological entities, confusion may arise between the two colors which show in the media.

Further, it is desirable to be able to identify and differentiate other closely related organisms, such as members of the families Aeromonaceae, Vibrionaceae, and Salmonella and Shigella spp. For example, the genus *Aeromonas* is closely related to coliforms and gives an almost identical biochemical test pattern. Further, the genus *Vibrio* is also an

5

important type of bacteria that grows under the same general conditions as coliforms. It is known to distinguish *Aeromonas* colonies from general coliforms by testing all colonies in a given sample for the presence of cytochrome oxidase. Undesirably, however, this requires another set of tests. Further, *Aeromonas* is common in water and food, and it is commonly indicated in test samples as general coliforms, which results in high a false positive error for general coliforms by most current test methods. The *Aeromonas* can be prevented from interfering with the coliform results by adding certain antibiotics to the medium. However, the antibiotic amounts added must be carefully controlled. Further, the antibiotics which have been conventionally used have short life spans in the media so that they lose their potency quickly in other than a frozen condition. It may often be desirable to be able to culture, identify and enumerate *Aeromonas* spp. which cannot be done if they are inhibited.

Further, in those cases where it is desirable to inhibit *Aeromonas*, it is desirable for a better method of so doing, one in which the shelf life of the medium is not appreciably reduced by the inclusion of an inhibitor.

Additionally, it is also desirable to distinguish strains of *Salmonella* and *Shigella* from *E. coli*, general coliforms and *Aeromonas*.

SUMMARY OF THE INVENTION

The present invention overcomes the disadvantages of prior art methods by providing a test method and medium for quantitatively or qualitatively identifying and differentiating biological entities in a test sample that may include a plurality of different biological entities.

The present invention introduces the use of "nonchromogenic" substrates to enhance the distinction among multiple colors produced by distinct biological entities present in the inventive test medium. Unexpectedly, it has been discovered that other "chromogenic" substrates present in the inventive test medium do not interfere with the substantially black color achieved with the nonchromogenic substrate. That is to say, so long as a given biological entity is responsive to the nonchromogenic substrate, aggregations thereof present in the test medium will show as a substantially black color—independent of whether such biological entity is responsive to one, two or more chromogenic substrates which are also present in the medium. The present invention exploits this hitherto unexplored property of nonchromogenic substrates.

In one form thereof, the present invention provides a test medium for detecting, identifying and qualifying or quantifying first and second biological entities. The test medium includes a nutrient base medium including ions of a salt, a chromogenic substrate and a nonchromogenic substrate. The first biological entity is responsive to the nonchromogenic substrate whereas the second biological entity is responsive to the chromogenic substrate. In this test medium, aggregations of the first biological entity present in the test medium are substantially black and aggregations of the second biological entity present in the test medium are a second color, the second color being distinguishable from the substantially black.

6

In a preferred form, the inventive test medium accounts for the first biological entity being responsive to the chromogenic substrate in addition to the nonchromogenic substrate. In such event, aggregations of the first biological entity present in the test medium will nonetheless show as substantially black.

Significantly, even though the aggregations of the first biological entity are responsive to both the first and second substrates in the preferred form, these aggregations still show as substantially black in the test medium. That is, the chromogenic substrate does not interfere with the substantially black color. Advantageously, this property of nonchromogenic substrates allows several different biological entities to be identified and differentiated in a single medium, aggregations of each biological entity having a visually distinguishable color.

In another preferred form of the above-described inventive medium, the medium further includes the antibiotic nalidixic acid to inhibit the growth of *Aeromonas*, spp. Advantageously, it has been found that nalidixic acid, as compared with cefsulodin, does not significantly reduce the shelf life of the test medium incorporating it.

In this connection, another form of the present invention provides a method of making a test medium for detecting at least one first type of biological entity and inhibiting a second type of biological entity from growing in the medium. The method includes the steps of combining desired substrates with a nutrient base medium; adding an inhibitor to the medium; and then sterilizing the medium by subjecting the medium to at least 100° C. Because the inhibitor is added as an initial step, subsequent sterile addition of inhibitor is unnecessary.

In another form thereof, the present invention provides a test medium for detecting, identifying and qualifying or quantifying first, second and third biological entities. The test medium includes a nutrient base medium including ions of a salt. First and second chromogenic substrates and a nonchromogenic substrate are provided in the test medium. The first and second biological entities are responsive to the first and the second chromogenic substrates, respectively, and the third biological entity is responsive to the nonchromogenic substrate. Aggregations of the first biological entity present in the test medium are a first color, aggregations of the second biological entity present in the test medium are a second color, and aggregations of the third biological entity present in the test medium are substantially black.

In a preferred form, the inventive test medium accounts for the third biological entity being responsive to the first and/or the second chromogenic substrates in addition to the nonchromogenic substrate. In such event, aggregations of the third biological entity will nonetheless show as substantially black.

It should be appreciated that the use of a nonchromogenic substrate along with one or more chromogenic substrates synergistically increases the number of biological entities that can be detected and distinguished in a single medium and synergistically increases the possible color combinations for a given set of biological entities to be detected. Stated another way, including a nonchromogenic component as one of the substrates synergistically increases the degrees

of freedom in selecting other substrates and corresponding colors for a test medium. This is so because an aggregation of the biological entity which is responsive to the nonchromogenic substrate will dependably show as substantially black. No combined color effects need be accounted for with the nonchromogenic substrates. For example, in a test medium including three chromogenic substrates and a nonchromogenic substrate, at least three combined color combination effects are avoided by using the one nonchromogenic component, as compared with using four chromogenic components.

The present invention, in another form thereof, provides a test medium capable of detecting, quantifying, and differentiating general coliforms and/or *E. coli* spp. under ambient light. The test medium comprises a nutrient based medium including a salt. A first substrate capable of forming a first water insoluble component of a first color in the presence of *E. coli* and the ions of the salt is provided in the medium. The first color is substantially black. A second substrate capable of forming a second water insoluble component of a second color in the presence of general coliforms is provided. The second color is visually distinguishable from the first color. Thus, colonies of *E. coli* present in the test medium are indicated by the first substantially black color and colonies of general coliforms are indicated by the second color.

In a preferred form of the above invention, the test medium further includes a third substrate capable of forming a third water insoluble component of a third color in the presence of *Salmonella* or *Shigella* spp. The third color is distinguishable from the first and second colors, whereby the test medium is capable of quantifying and/or differentiating *E. coli*, general coliforms and *Salmonella* or *Shigella* spp. Further, the substrates are selected such that general coliforms present in the test medium will also react with the third substrate to form a water insoluble component which includes the third color. Consequently, general coliform colonies are indicated in the test medium as a fourth color, the fourth color being a combination of the second color and the third color. The fourth color is visually distinguishable from the first and third colors. Still further, the substrates can be selected such that *Aeromonas* spp. form an insoluble component of the second color by reacting with the second substrate, but not the first and third substrates. Thus, in the inventive test medium, *E. coli* colonies will be generally black, general coliform colonies will be the fourth color, *Aeromonas* colonies will be the second color and *Shigella* or *Salmonella* colonies will be the third color.

In another form thereof, the present invention provides a method for detecting, quantifying and differentiating under ambient light general coliforms, *E. coli*, and at least one of the genera *Aeromonas*, *Salmonella* or *Shigella* in a test sample. The method comprises the steps of providing a nutrient base medium including first, second and third substrates. Each of the substrates is capable of forming a water insoluble component in the presence of at least one of general coliforms, *E. coli*, *Aeromonas*, *Salmonella* or *Shigella*. The substrates are selected such that colonies of *E. coli* produced in the test medium are a first color, colonies of general coliforms produced in the test medium are a second color, and colonies of one of *Aeromonas* and *Salmonella* or

Shigella produced in the test medium are a third color. Each of the colors are visually distinguishable. The test medium is inoculated with the test sample and then incubated. The test medium is then examined for the presence of first colonies having the first color, second colonies having the second color, and third colonies having the third color. The first colonies are *E. coli*, the second colonies are general coliforms, and the third colonies are one of *Aeromonas*, *Salmonella* or *Shigella*.

In a preferred form thereof, the inventive method further includes adding ions from a salt to the test medium to react with one or more of the substrates. In so doing, a precipitate is produced which shows as a substantially black color in the presence of the specific enzyme for that substrate. A preferred compound for forming the substantially black color in the presence of the ions of the salt consists of a β -D-glucuronide. These compounds release an aglycone when hydrolyzed which forms a substantially black insoluble complex in the presence of ions.

In another preferred form of the inventive method, the method further comprises examining the test medium for the presence of fourth colonies having a fourth color, wherein the substrates are selected such that colonies of *Aeromonas* are the third color and colonies of *Salmonella* or *Shigella* are the fourth color, the fourth color being visually distinguishable from the first, the second and the third colors. More preferably, the substrates are selected such that the first color is substantially black, the second color is substantially blue-violet, the third color is substantially red-pink and the fourth color is substantially teal-green.

In another preferred form of the inventive method, the substrates are selected such that colonies of *Aeromonas* as well as colonies of *Plesiomonas* and *Vibrios* are indicated as the third color.

One advantage of the present invention is that it uses a nonchromogenic substrate along with one or more chromogenic substrates and thereby synergistically increases the degrees of design freedom in selecting colors for the inventive test medium. This is so because the chromogenic substrates do not interfere with the substantially black precipitate formed by the nonchromogenic substrate.

Another advantage of the present invention is that enables the quantification, identification and differentiation of four (4) different bacterial strains simultaneously in a single test medium using a single test sample, under ambient lighting. Subsequent tests with their concomitant extra time spent and extra costs are avoided. Of course, the inventive test medium of the present invention could also be used purely for qualitative purposes, as a mere presence/absence (P/A) test.

Yet another advantage of the present invention is that the substrates are selected such that the colors are easy to visually distinguish from one another without the need for UV light or other visual aids, other than, perhaps, magnification means. For example, in a preferred embodiment, *E. coli* colonies are clearly indicated by a precipitate having a substantially black color, general coliform colonies are indicated by a blue-violet color, *Aeromonas* colonies are indicated by a red-pink color, and *Salmonella* or *Shigella* colonies are indicated by a teal-green color. Because these colors are visually so distinct, confusion among the colors is greatly reduced as compared to prior art media.

Another advantage of the test medium of the present invention is its flexibility and ease of use. The incubation temperature is not critical as growth and differentiation of the organisms mentioned may occur within an optimum range. Therefore, resuscitation steps are avoided and inhibition of temperature sensitive strains does not occur. Also, inexpensive equipment may be used.

Yet another advantage of the present invention is that it intensifies the color distinction obtained in a test medium for identifying and differentiating *E. coli* from general coliforms. In a preferred test medium, *E. coli* colonies present a substantially black color, whereas general coliforms present a red-pink color, the distinction therebetween being much more apparent than in prior art test media. Confusion between the two colors is therefore greatly reduced.

Still another advantage of the present invention is that it enables the identification and differentiation of *Aeromonas* spp. from general coliforms. Prior art test media undesirably require using a cefsulodin inhibitor for preventing *Aeromonas* spp. from growing therein. However, the use of cefsulodin as an inhibitor requires an extra step in the process, viz., sterile addition of filter sterilized antibiotic, and is difficult to control. Further, the presence of cefsulodin significantly reduces the effective shelf life of the medium. Further, the use of an inhibitor, obviously, prevents the detection and quantification of *Aeromonas* spp. Advantageously, with the present invention, *Aeromonas* spp. can be detected, quantified and differentiated from general coliforms in a single medium.

As a related advantage, if it is nonetheless desired to inhibit colonies of *Aeromonas* spp. from growing in the test medium, the present invention provides a superior means for doing so. Specifically, preferred forms of the present invention employ nalidixic acid as an inhibitor, which has been shown to have a far less deleterious effect in the shelf-life of the medium incorporating it. Further, nalidixic acid can be added as part of the initial medium formulation prior to sterilization, thereby avoiding a costly and difficult process step which is required with cefsulodin. Finally, nalidixic acid is much less expensive than cefsulodin.

Another advantage of the present invention is that it can provide a test medium for qualitative or quantitative testing. That is, the test media in accordance with the present invention can be used as mere presence/absence test devices, or can be used to quantify various biological entities showing as different colored colonies on the inventive test media.

DETAILED DESCRIPTION OF THE INVENTION

The method and medium of the present invention allow the simultaneous detection, quantification, identification and differentiation of a variety of selected biological entities in a sample of mixed populations of biological entities. The inventive method and medium are particularly useful for the detection, quantification, identification and differentiation of *E. coli* and general coliforms, and further quantitative identification and differentiation of other selected biological entities, including *Aeromonas*, *Salmonella*, *Shigella*, *Pseudomonas*, and *Vibrio* bacterial species.

The method and test media incorporating the present invention utilize the fact that the enzymatic activity of

biological entities and specifically of bacteria varies with the genus, and/or family of bacteria of interest. The method and test media incorporating the present invention further utilize the fact that various enzyme identifying substrate complexes can be used to identify specific enzymes with the production of distinctive colors. Significantly, the method and test media incorporating the present invention exploit the fact that chromogenic substrates present in a test medium do not interfere with the substantially black color produced by nonchromogenic substrates.

While nonchromogenic substrates are known in the art, per se, their distinct properties vis-à-vis chromogenic substrates have been unrecognized. However, the behavior of a nonchromogenic substrate in a medium including a combination of chromogenic substrates is unique. To illustrate, aggregations of a biological entity which is responsive to two chromogenic substrates will typically show in a test medium as a combination of the two colors produced upon cleavage of the two respective substrates. If three chromogenic substrates are involved, the combined color effect could be prohibitively difficult to predict and account for. Further, inherent variations in the amount of enzymes produced by particular strains of biological entities can result in different shades or hues of colors upon cleavage of the chromogenic substrates. Consequently, the colors can be difficult to distinguish for the lay person examining the test medium. Chromogenic substrates must therefore be chosen in view of the other chromogenic substrates planned for inclusion in a given test medium.

Such is not the case with the nonchromogenic components. While aggregations of biological entities which are responsive to chromogenic substrates in addition to nonchromogenic substrates may show in the test medium as having a colored or fluorescent "halo," such aggregations nonetheless appear substantially black and are therefore easy to identify. Unlike chromogenic substrates, multiple "degrees of freedom" are achieved with the nonchromogenic components by not having to take into account combined color effects.

Using a nonchromogenic substrate enables a single test medium to differentiate four (4) different bacterial strains with four (4) visually distinguishable colors. The black color is superior in that it is difficult to mistake. Further, the substantially black pigmentation does not diffuse so that the location of the colonies is precisely known and the colonies can be accurately counted. The nonchromogenic substrates produce an insoluble chelated compound which is different than the dimer which is produced by the chromogenic substrates.

The inventive test medium and method allows not only a detection, quantification or qualitative identification and differentiation of general coliforms and *E. coli*, but also of *Salmonella*, *Shigella* and *Aeromonas*, as well as *Plesiomonas* and *Vibrio*. *Plesiomonas* and *Vibrios* species are determined but not differentiated from *Aeromonas* species as they are very closely related.

DEFINITIONS

Biological entities, such as general coliforms, *E. coli*, etc., are herein referred to as being "responsive" to certain

chromogenic and nonchromogenic substrates. More specifically, a biological entity will predictably produce specific enzymes when the entity is present in a test medium such as the one described hereinbelow. These enzymes will selectively cleave chromogenic and nonchromogenic substrates. Upon cleavage, these substrates produce a color in the test medium. The mechanism for producing the color is different for chromogenic and nonchromogenic substrates, as described hereinbelow.

Microorganisms having β -galactosidase activity include those commonly known by the designation "coliform." There are various definitions of "coliform," but the generally accepted ones include bacteria which are members of the Enterobacteriaceae family, and have the ability to ferment the sugar lactose with the evolution of gas and acid. Most coliforms are positive for both α - and β -galactosidase. That is, they produce both α - and β -galactosidases.

Microorganisms having β -glucuronidase activity in addition to galactosidase activity primarily include most strains of *Escherichia coli*. That is, *E. coli* is positive for both α - and β -galactosidase as well as β -glucuronidase.

The term "general coliforms" as used in this application refers to coliforms other than the various strains of *E. coli*. These "general coliforms" are gram-negative, non-spore forming microorganisms generally having α - and β -galactosidase activity (i.e., lactose fermenters), but not having β -glucuronidase activity, and having the ability to ferment the sugar sorbitol.

For purposes of this specification, a "chromogenic substrate" is a substrate which needs no additional chemicals present in the test medium upon hydrolysis for color production. That is, a chromogenic substrate is cleaved by the specific enzyme corresponding to that substrate to form a dimer with the color being concentrated in the area of cleavage of the substrate. Many chromogenic substrates are known in the art. For purposes of this specification "chromogenic" includes fluorogenic substrates. The products of fluorogenic substrates require ultraviolet (UV) light to be detected and are more soluble than preferred chromogenic substrates, so are therefore generally not preferred for use with the test media disclosed hereinbelow.

Certain substrates, referred to herein as "nonchromogenic," produce a dark, substantially black precipitate in the presence of ions of a salt and enzyme activity. For example, 8-hydroxyquinoline- β -D-glucuronide, when included in a medium along with a salt that produces ions, such as ferric ammonium citrate, will produce a substantially black precipitate in the presence of β -glucuronidase produced by *E. coli* or other biological entities. More specifically, upon cleavage of the nonchromogenic substrate by the particular enzyme, a substantially black water-insoluble complex forms in the medium. The substantially black precipitate consists of the ferric ions and the aglycone released when the substrate is hydrolyzed by the glucuronidase from *E. coli*. This precipitate is a chelated compound which does not diffuse. Nor is the substantially black color susceptible to interference from chromogenic compounds present in the test medium.

For purposes of this specification a "nonchromogenic substrate" means that a chemical in addition to those used

with chromogenic components must be present in the test medium when the substrate is cleaved by its corresponding enzyme. The substantially black precipitate formed thereby is a combination of the substrate-salt complex and is not a dimer as is formed by the "chromogenic compounds."

For purposes of this specification, the expression "under ambient light" refers to the visible spectrum, i.e., colors which can be seen and distinguished with the naked eye. A colony present in a test medium which requires ultraviolet light to be seen, for example, would not fall under the definition "under ambient light". However, it is to be understood that the term "under ambient light" includes using a magnification device, if necessary. Magnification can be especially helpful when counting numerous colonies. The term "visually distinguishable" refers to two or more colors which can be differentiated under ambient light.

For purposes of this specification, the term "substantially black" includes dark brown to black, and also includes black with various colored halos, such as red-violet, green, fluorescent, etc.

For further purposes of this specification, color names recited herein are given as guidance, but it is to be understood that the color names are to be read broadly. That is, there can be overlap among the recited colors. This is because, as discussed, biological entities produce varying amounts of enzymes, which in turn affects the shade or hue of the resulting color.

The term " β -galactosidase substrate" as used herein refers to a β -galactoside comprising galactose joined by β -linkage to a substituent that forms a water insoluble colored compound when liberated by the action of β -galactosidase on the substrate. Similarly, the term " α -galactosidase substrate" as used herein refers to α -galactoside comprising galactose joined by α -linkage to a substituent that forms a water insoluble colored compound when liberated by the action of α -galactosidase on the substrate. The term " β -glucuronidase substrate" as used herein refers to a β -glucuronide comprising glucuronic acid joined by β -linkage to a substituent that forms a water insoluble colored precipitate when liberated by the action of β -glucuronidase on the substrate.

The α - and β -galactosidase substrates and compounds and any other substrates described herein as well as the β -glucuronidase substrates and compounds and any other substrates described herein may comprise carboxylate salts formed by reacting a suitable base with the appropriate galactoside or glucuronic carboxyl group. Suitable bases include alkali metal or alkaline earth metal hydroxides or carbonates, for example, sodium hydroxide, potassium hydroxide, calcium hydroxide, magnesium hydroxide, and corresponding carbonates; and nitrogen bases such as ammonia, and alkylamines such as trimethylamine, triethylamine and cyclohexylamine.

Designing a Test Medium for Specific Biological Entities

Certain members of the family Enterobacteriaceae can be distinguished by the presence of α -galactosidase activity in the absence of β -galactosidase activity, or vice-versa. For example, most *Salmonella* and *Shigella* spp. are positive for α -galactosidase, but negative for β -galactosidase. Similarly,

Aeromonas spp. can be distinguished from other members of the family Enterobacteriaceae by the presence of β -galactosidase activity in the absence of α -galactosidase activity. The method and medium incorporating the present invention are designed to take advantage of these distinguishing characteristics. For example, the specificity of enzyme activity for *Salmonella* and *Aeromonas* spp., as opposed to general coliforms, can be exploited, as illustrated below.

The method described herein is particularly suitable for the detection, quantification or qualitative identification and differentiation of the different classes of microorganisms described previously, viz., general coliforms, *E. coli*, *Aeromonas* and *Salmonella* and *Shigella* spp. Although the inventive method is particularly suitable for these particular microorganisms, it is not limited thereto. Instead, the techniques described herein have application to the identification and differentiation of a wide variety of biological entities.

That is, specific biological entities are "responsive" to various substrates. More particularly, these biological entities predictably produce or contain known enzymes. Substrates, either chromogenic or nonchromogenic, can be selected which, in the presence of a particular enzyme(s), will form an insoluble component of a predictable and distinguishable color. Multiple substrates can be selected to simultaneously identify a plurality of distinct biological entities in a single test medium, aggregations of each distinct entity being identifiable by a separate, distinguishable color. Further, while the preferred embodiments disclosed herein distinguish all of the various aggregations present in a test medium under ambient light, as that term is defined herein, such is not necessary. For example, several substrates disclosed herein require the use of ultraviolet light for the aggregations present in the medium to be seen.

Table I lists various enzymes whose presence may be detected using certain of the substrates listed in Table II.

TABLE I

Enzymes and Abbreviations	
Aara = α -D-arabinopyranosidase	Bglu = β -D-glucopyranosidase
Agal = α -D-galactopyranosidase	Bgluc = β -D-glucuronidase
Aglu = α -D-glucopyranosidase	Bmaa = β -D-mannopyranosidase
Bcel = β -D-cellobiosidase	Bxyl = β -D-xylopyranosidase
Rfuc = β -D-fucopyranosidase	Nagal = N-acetyl- β -D-galactopyranosidase
Bgal = β -D-galactopyranosidase	Naglu = N-acetyl- β -D-glucopyranosidase
Afuc = α -D-fucopyranosidase	Bara = β -D-arabinopyranosidase
Bxyl = β -D-xylopyranosidase	Acel = α -D-cellobiosidase
Aman = α -D-mannopyranosidase	Agduc = α -D-glucuronidase
Axyl = α -D-xylopyranosidase	Nagduc = N-acetyl- β -D-glucuronidase
esterase	

TABLE II

Various Substrates and Color Upon Cleavage	
6-chloro-3-indolyl substrates	Pink
5-bromo-4-chloro-3-indolyl substrates	Teal
3-indolyl substrates	Teal
N-methylindolyl substrates	Green
nitrophenyl substrates	Yellow
nitroaniline substrates	Yellow

TABLE II-continued

Various Substrates and Color Upon Cleavage	
8-hydroxyquinoline substrates (and ion of salt)	Substantially black
cyclohexenoescluletin substrates (and ion of salt)	Substantially black
esculetin substrates (and ion of salt)	Substantially black
quinoline substrates (and ion of salt)	Substantially black
5-Iodo-3-Indolyl substrates	Purple
5-Bromo-6-Chloro-3-Indolyl substrates	Magenta
6-Fluoro-3-Indolyl substrates	Pink
coumarin substrates	Fluorescent
fluorescein substrates	Fluorescent
rhodamine substrates	Fluorescent
resorufin substrates	Fluorescent

Specific substrate compounds applicable for use with the test medium of the present invention are available as follows:

5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside (X-gal) is a commercially available β -galactosidase substrate that produces an insoluble precipitate having an approximately teal color when reacted upon by β -galactosidase and is available from Biosynth International, Naperville, Ill.

6-chloro-3-indolyl- β -D-glucuronide is a compound which produces an insoluble precipitate having a magenta color, the preparation of which is described in the aforementioned incorporated by reference U.S. Pat. No. 5,210,022 and is available from Research Organics, Cleveland, Ohio.

The compound 5-bromo-4-chloro-3-indolyl- β -D-glucuronide (X-gluc) is a commercially available β -glucuronide that produces an insoluble precipitate having an approximately teal color when reacted upon by β -glucuronidase. Similarly, indoxyl- β -glucuronide is a similar compound, the preparation of which is described in the aforementioned article by Ley et al., in *Can J Microbiol.*, the disclosure of which is incorporated by reference.

Another suitable β -galactoside is the compound 6-chloro-3-indolyl- β -D-galactoside which produces an insoluble precipitate having a pink/magenta color, the preparation of which is described in the aforementioned U.S. Pat. No. 5,210,022.

Other suitable compounds applicable as substrates in the practice of the present invention are specified in U.S. Pat. No. 5,210,022, all of which are incorporated herein by reference.

The substrate 8-hydroxyquinoline- β -D-glucuronide is a commercially available β -glucuronide that, in the presence of metallic ions such as iron, produces an insoluble precipitate having a substantially black color when reacted upon by β -glucuronidase and in the presence of other α - or β -galactoside substrates. 8-hydroxyquinoline- β -D-glucuronide is available from Biosynth International, Naperville, Ill.

Further, a salt providing ions suitable for use with the present invention is ferric ammonium citrate, available from Sigma Chemical, St. Louis, Mo. The cyclohexenoescluletin substrates are described in James et al., *Appl. & Envir. Micro.* 62:3868-3870 (1996) and in the presence of ferric ions, produce an insoluble substantially black precipitate.

N-methyl-indolyl substrates such as N-methylhydroxy- β -D-galactopyranoside are commercially available from Biosynth International, Naperville, Ill.

Nitrophenyl substrates, such as 2-nitrophenyl- β -D-galactopyranoside, are commercially available from Biosynth International, Naperville, Ill. Similarly, nitroaniline compounds are available for synthesis through Sigma Chemical, St. Louis, Mo.

Other substrates producing a substantially black color include esculetin substrates such as cyclohexenoesculetin- β -D-galactoside, which is described in James et al., *Appl. & Envir. Microbiol.* 62:3868-3870 (1996). Quinoline substrates, such as 8-hydroxyquinoline- β -D-galactopyranoside and 8-hydroxyquinoline- β -D-glucuronide are available through Biosynth International, Naperville, Ill.

Iodo-indolyl substrates, such as 5-iodo-3-indolyl- β -D-galactopyranoside are available through Biosynth International, Naperville, Ill.

Several fluorescent substrates are suitable for use with the present invention. Coumarin substrates such as 4-methylumbelliferyl substrates and 5-trifluoromethylumbelliferyl substrates are commercially available from Biosynth International, Naperville, Ill. Also suitable are fluorescein substrates, rhodamine substrates, and resorufin substrates. No commercial source is known for these three substrates but components are available from Sigma Chemical, St. Louis, Mo.

While specific examples of substrates suitable for use with the present invention have been enumerated hereinabove, such is not to be construed as limiting the invention in any manner. Instead, one of ordinary skill in the art can use Table IV and V hereinbelow to identify a virtually limitless number of substrates.

Preparation of Test Medium

The test medium is formed by combining the desired substrates with a nutrient base medium. The nutrient base medium can be any culture medium known in the art for providing the maintenance and reproduction of living cells. Generally, such media include nutrients, buffers, water, and sometimes a gelling agent. Possible gelling agents include agars, pectins, carrageenans, alginates, locust bean, and xanthins, among others.

The following is an example of the preparation of a test medium suitable for use in this invention. This example coincides with Example I, below.

The substrates 8-hydroxyquinoline- β -D-glucuronide, 5-Bromo-4-chloro-3-indolyl- α -D-galactopyranoside, and 6-Chloro-3-indolyl- β -D-galactopyranoside are added in quantities of 250 mg/L medium; 70 mg/L medium; and 175 mg/L medium, respectively. The substrates are added directly to the hot (75°-85° C.) medium (formula below) in a blender prior to sterilization.

Standard agar medium may be made by adding 15 gm of bacteriological quality agar gum to the following nutrient formula

Pancreatic Digest of Casein	5.0 gm
Yeast Extract	3.0 gm

-continued

Dipotassium Phosphate	.3 gm
Deionized Water	990 ml
Ferric Ammonium Citrate (sterilized separately from the other components)	800 mg in 10 ml deionized water

and then sterilizing at 121° C. for 15 minutes. The medium should be adjusted to result in a pH of 7.0. The sterilized agar medium is allowed to drop to a temperature of 45° C. in a water bath and then the sterile solution containing the substrates prepared as described above is added. The medium is mixed thoroughly and poured into sterile petri plates at a volume of 20 ml/plate.

A pectin-based test medium may be prepared using the same steps described above except that 25 gm of low methoxyl pectin is used as the solidifying agent and the medium is poured at room temperature into petri plates containing a thin gel layer containing calcium ions which combine with the pectin to form a solid gel. A suitable pectin culture medium is described in U.S. Pat. No. 4,241,186 and U.S. Pat. No. 4,282,317, the disclosures of which are incorporated herein by reference. A pectin-based medium is preferred over a standard agar medium because it has the advantages of convenience and temperature independence for the user. The use of pectin media is well described and accepted as a result of AOAC collaborative studies and other published and in-house investigations.

A suitable pectin medium is commercially available from Micrology Laboratories, LLC under the trademark Easy-gel®. Aqueous based medium without gelling agent is available from Micrology Labs, Gosben Ind., for use with membrane filters.

Inoculation of the Test Medium With the Sample

The test medium may be inoculated with a sample to be tested for the presence of microorganisms by any method known in the art for inoculating a medium with a sample containing microorganisms. For example, the sample to be tested may be added to the petri plates prior to adding the nutrient base medium (pour plate technique) or spread on the surface of the plates after they have cooled and solidified (swab or streak plate technique). Liquid samples may also be filtered through a micropore (0.45 micrometer size) membrane filter which is then placed on the surface of a solid medium or on a pad saturated with the medium.

Incubation of the Test Medium

The inoculated test medium is incubated for a sufficient time and at such a temperature for individual microorganisms present in the sample to grow into detectable colonies. Suitable incubation conditions for growing microorganisms in a medium are known in the art. Commonly, the test medium is incubated for about 24-48 hours at a temperature of about 30°-40° C.

Unless inhibitors of the general microbial population are used, the general microbial population as well as general coliforms, *E. coli*, *Aeromonas* spp., and *Salmonella* spp. and *Shigella* spp. will grow in the incubated test medium. Because the precipitates formed are insoluble in the test medium, they remain in the immediate vicinity of micro-

ganisms producing the various enzymes. As the microorganisms reproduce to form colonies, the colonies show as colony forming units having the color produced by the particular substrate.

For example, *E. coli* produces β -galactosidase and α -galactosidase, but, unlike general coliforms and *Aeromonas* spp., also produces β -glucuronidase. Therefore, insoluble precipitates of each of the β -galactosidase substrate, the α -galactosidase substrate and the β -glucuronide substrate are formed by the action of the respective enzymes such that colonies of *E. coli* show as a substantially black color, sometimes having a violet-blue halo therearound.

General coliforms produce β -galactosidase and α -galactosidase and consequently cleave both the α -galactosidase and β -galactosidase substrates. In the present example, the 5-Bromo-4-chloro-3-indolyl- α -D-galactoside substrate produces a blue-green or teal color, whereas the 6-Chloro-3-indolyl- β -D-galactoside produces a pink, or red-pink color. Thus, general coliform colonies will show as a blue-violet color, which is a combination of the colors produced by each of the α - and β -galactosides, respectively.

Significantly, however, it has been found that *Aeromonas* spp., which are closely related to coliforms, and give an almost identical biochemical test pattern, are β -galactosidase positive and α -galactosidase negative. That is, *Aeromonas* spp. will not hydrolyze the α -galactoside substrate. Therefore, *Aeromonas* colonies present in the test medium will show as colonies having the pink-red color produced by the β -galactoside substrate.

Further, it has been found that most members of the genera *Salmonella* and *Shigella* are positive for α -galactosidase, but negative for β -galactosidase. That is, *Salmonella* and *Shigella* will not hydrolyze the β -galactosidase substrate. Therefore, colonies of *Salmonella* and *Shigella* present in the test medium will appear as a teal, or blue-green color produced by the α -galactoside substrate. Occasionally, *Shigella* colonies will appear black with a blue-green halo since some strains of *Shigella* are positive for β -glucuronidase, and some strains of *Shigella* will appear blue/purple since some unusual strains are positive for both α -galactosidase and β -galactosidase.

Examination of the Test Medium and Enumeration of Microorganisms

The substrates selected for the above example produce three distinct colors, and general coliforms are indicated by a fourth color which is a combination of two of the three colors. That is, *E. coli* colonies show as substantially black, general coliform colonies show as blue-violet, *Aeromonas* colonies show as red-pink, and *Salmonella* and *Shigella* colonies show as teal-green. While the individual shades of these colors may vary somewhat in the test medium due to factors such as varying enzyme production of the biological entities, it has been found that these four colors are distinct enough so that confusion amongst them is unlikely.

The colonies of each type of microorganism may be enumerated by counting the colonies or by other methods known in the art for enumerating microorganisms on a test

plate. The number of colonies of each type indicates the number of microorganisms of each type originally present in the sample before incubation.

Optional Ingredients

Inhibitors

The method of the present invention does not require inhibitors. However, the medium may be made more selective for general coliforms and *E. coli* if desired by the addition of various compounds that are inhibitory to the general microbial population, but have little or no effect on coliforms. Following are some compounds which may be used: a) bile salts, about 0.3 g/liter, b) sodium lauryl sulfate, about 0.2 g/liter, c) sodium desoxycholate, about 0.2 g/liter, d) Tergitol 7, about 0.1 ml/liter. The use of one or more of these compounds reduces the background (non-coliform) microorganism presence and makes a less cluttered plate and eliminates the possibility of inhibition or interference by the non-coliform organisms in the sample. The use of certain antibiotics may accomplish the same result.

Cefsulodin is commonly used in currently available test media to inhibit *Aeromonas* spp. However, the use of cefsulodin as an inhibitor requires an extra step in the process, viz., sterile addition of filter sterilized antibiotic. This step is difficult to control. Further, the presence of cefsulodin significantly reduces the effective shelf life of the medium. It has been found that Nalidixic acid can be used instead of Cefsulodin to inhibit *Aeromonas* spp. with about the same efficacy. Nalidixic acid is preferable because it can survive the approximately 120° C. temperature reached in autoclaving the test media. Therefore, unlike cefsulodin, nalidixic acid can be added to the test media as part of the initial media formulation prior to sterilization (see, preparation of test medium, above). It also follows that the resistance of the nalidixic acid to unfavorable environmental conditions will result in a longer shelf life for a medium containing it as compared to cefsulodin.

Inducers

It is possible that the enzyme production of the general coliforms may be enhanced by the addition to the medium formulations of very small amounts of substances known as enzyme inducers. One specific inducer for β -galactosidase is available and is known chemically as isopropyl- β -thiogalactopyranoside. Adding approximately 100 mg/liter of medium has a positive and noticeable effect on the speed of enzyme production for some species of coliforms. Other enzyme inducers are available and may be added to media formulations if enhanced enzyme production is deemed helpful.

EXAMPLES

Listed below are broad examples of test media enzyme substrate combinations to be used in combination with the nutrient formula discussed above or other suitable nutrient formulas which may be prepared in practicing the present invention.

Table III illustrates the flexibility of the preferred embodiments incorporating the present invention. Table III is a matrix of some of the possible four-color combinations

available for the preferred biological entities *E. coli*, general coliforms, and at least one of the genera *Aeromonas*, *Salmonella* or *Shigella* to be detected by using the teachings of this disclosure. Other color combinations are possible. In many cases, a plurality of different substrates will achieve a desired result, the only difference being the colors detected for a specific enzyme. The preferred color choice for the detection of *E. coli* is denoted with an asterisk in Table III, depending on the colors chosen to detect other microorganisms. As discussed above, the substantially black color is preferred because other chromogenic substrates do not interfere with it and the substantially black color is easy to distinguish from the other colors.

As discussed above, the use of Table III requires taking into account the combined color effect discussed above which is produced by the inclusion of multiple chromogenic substrates in a single medium. For example, with reference to the first entry in Table III, it can be understood that general coliforms will appear as a combination of (1) red-pink (magenta) and (2) teal, the resulting color being blue-violet. This is the case because general coliforms are responsive to two chromogenic substrates. Similarly, general coliforms will show in a test medium in accordance with the second entry of Table III as a combination of (1) red-pink (magenta) and (2) yellow.

TABLE III

Color possibilities for detection of preferred microorganisms								
desired color	red-pink or magenta	teal	green	yellow	black	fluorescent	fluorescent	fluorescent
1	general coliforms <i>Aeromonas</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i> *	<i>E. coli</i>		
2	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i> *	<i>E. coli</i>		
3	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i> *	<i>E. coli</i>		
4	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i> *	<i>E. coli</i>		
5	<i>E. coli</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i> *	<i>E. coli</i>		
6	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i> *	<i>E. coli</i>		
7	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	
8	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i> *	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	
9	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i> *	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	
10	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i> *	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	
11	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i> *	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>	
12	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>		
13	<i>F. coli</i>	general coliforms <i>Aeromonas</i>	<i>F. coli</i>	<i>F. coli</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>F. coli</i>		
14	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	<i>E. coli</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>		
15	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	general coliforms <i>Salmonella</i> / <i>Shigella</i>	<i>E. coli</i>		
16	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i>	<i>E. coli</i> *	<i>E. coli</i>	general coliforms <i>Aeromonas</i>	general coliforms <i>Salmonella</i>

TABLE III-continued

Color possibilities for detection of preferred microorganisms							
desired color*	red-pink or magenta	teal	green	yellow	black	fluorescent	fluorescent
						monas	/shigella

* = preferred color for *E. coli*

Table IV is a partial list of enzyme patterns for biological entities preferred to be detected in accordance with the teachings of this disclosure. It is to be understood that one of ordinary skill in the art would readily recognize that other enzymes which are known and have been produced, and enzymes which are known only on a theoretically level, would also perform satisfactorily.

understood that Table V teaches a large quantity of substrates possible for use in accordance with the present invention. Many of the substrates identified by the above-described use of table V are commercially available, whereas the method for producing other identified substrates is described in the literature. Still other substrates identified by using table V are only theoretically possible.

TABLE IV

ENZYME NAME	<i>E. coli</i>	GENERAL COLIFORM	<i>Aeromonas</i>	<i>Salmonella/Shigella</i>	<i>Plesiomonas</i>	<i>Vibrio</i>
Aara=	+	+	—	+		
α -D-arabinopyranosidase						
Agal=	+	+	—	+		
α -D-galactopyranosidase						
Aglu=	—	+	+	—	+	
α -D-glucopyranosidase						
Bcel=	—	+	—	—	—	—
β -D-cellobiosidase						
Bfuc=	+	+	+	—	—	—
β -D-fucopyranosidase						
Bgal=	+	+	+	—	+	+
β -D-galactopyranosidase						
Bgal=	+	+	+	—	—	—
β -D-glucopyranosidase						
Bgluc=	+	—	—	—	—	—
β -D-glucuronidase						
Bman=	+	+	—	+	+	+
β -D-mannopyranosidase						
Bxyl=	—	+	—	—		
β -D-xyloripranosidase						
Nagal=	—	+	+	—	+	+
N-acetyl- β -D-galactopyranosidase						
Naglu = N-acetyl- β -D-glucopyranosidase	—	+	+	—	+	+
Aman=	—	—	—	—	—	—
α -D-mannopyranosidase						
esterase =	—	—	—	+	—	—
esterase						

Table V is a matrix which teaches a wide variety of substrates and their associated colors for use in test media in accordance with the teachings of this disclosure. The left hand side of Table V indicates the color that will result when the listed chromogenic component is cleaved from its corresponding substrate by the specific enzyme for that substrate. In the case of the nonchromogenic components, the color is substantially black and the reaction mechanism requires the presence of ions of a salt upon cleavage of the substrate, as explained above.

Test enzymes which are produced by certain biological entities (see Table IV) are found at the right hand side of table V. "Substrate components" are shown to the left of the specific test enzymes. Each of the substrate components listed on the right hand side of table V can be combined with any of the chromogenic or nonchromogenic components listed on the left hand side of table V to identify a specific substrate for use in a test medium. It can therefore be

Nonchromogenic components are included at the bottom left hand side of Table V, and are different from the chromogenic components because they do not form specific colors upon cleavage. Instead, the quinoline or esculetin components combine with ions of a salt (e.g., ferric salt) which must be present in the medium when the substrate is cleaved by the specific enzyme. The substantially black precipitate formed by the nonchromogenic components is a combination of the quinoline or esculetin—iron complex rather than a dimer which is formed by the chromogenic components.

Unlike nonchromogenic components, the chromogenic components should be selected in view of all other chromogenic components selected for the medium and in view of the enzyme patterns of the entities to be detected. The selection of chromogenic components should maximize the distinction among the respective colors produced.

While many various chromogenic component and substrate component/enzyme possibilities are taught by Table V, other possibilities within the scope of the appended claims would be possible by one of ordinary skill in the art. For example, as shown in Table V, one of ordinary skill in the art could combine an N-acetyl group with many of the sugars of the substrate components listed in Table V. For example, an N-acetyl group could be combined with β -D-mannopyranoside to form N-acetyl- β -D-mannosaminide, the corresponding enzyme being N-acetyl- β -D-mannosaminidase. Any of the chromogenic components or nonchromogenic components listed on the left hand side of Table V could then be combined with the substrate component to identify a substrate. If the substrate is commercially available or the method of making it is known, the substrate could be used in a test medium. Upon cleavage of the substrate by the corresponding enzyme in the test medium, the color listed will appear.

Generally, the teachings of this disclosure can be used as follows to make a test medium for detecting various microorganisms or cell types. First, the microorganisms desired to be detected and differentiated are selected. The preferred organisms to be detected are *E. coli*, general coliforms, and at least one of the genera *Aeromonas*, *Salmonella* or *Shigella*. Enzymes produced by the selected organisms can be identified with reference to Table IV. Equipped with knowledge of specific enzymes produced by each microorganism, one can then identify corresponding substrate components from the right hand side of Table V. Depending upon the color desired, one can select a chromogenic or nonchromogenic component from Table V to be combined with the substrate component to identify a substrate for inclusion in the test medium. If the substrate thereby identified is commercially available or the method of its synthesis is known, the substrate can be used in the test medium.

TABLE V

COLOR COMPONENT AND SUBSTRATE COMPONENT MATRIX		
CHROMOGENIC COMPONENT &	(COLOR)	SUBSTRATE COMPONENT -- TEST ENZYME
6-fluoro-3-indolyl-	(pink)	α -D-arabinopyranoside -- Aara.
6-chloro-3-indolyl-	(pink/red)	α -D-cellopyranoside -- Acel.
5-bromo-6-chloro-3-indolyl-	(magenta)	α -D-fucopyranoside -- Afuc.
3-indolyl-	(teal)	α -D-galactopyranoside -- Agal.
5-bromo-4-chloro-3-indolyl-	(teal)	α -D-glucuronide -- Agluc.
5-iodo-3-indolyl-	(purple)	α -D-mannopyranoside -- Aman.
N-methylindolyl-	(green)	α -D-xylopyranoside -- Axyl.
4-methylumbelliferyl-	(fluorescent)	β -D-arabinopyranoside -- Bara.
rhodamine-	(fluorescent)	β -D-cellopyranoside -- Bcel.
fluorescein-	(fluorescent)	β -D-fucopyranoside -- Bfuc.
resorufin-	(fluorescent)	β -D-galactopyranoside -- Bgal.
coumarin	(fluorescent)	β -D-glucopyranoside -- Bglu.
nitrophenyl-	(yellow)	β -D-glucuronide -- Bgluc.
nitroaniline	(yellow)	β -D-mannopyranoside -- Bman.
NONCHROMOGENIC COMPONENT	(COLOR)	β -D-xylopyranoside -- Bxyl.
		N-acetyl- β -D-galactosaminide -- Nagal
8-hydroxyquinoline plus ions-	(substantially black)	N-acetyl- β -D-glucosaminide -- Naglu
3,4-cyclohexenesculetin plus ions	(substantially black)	N-acetyl- β -D-glucuronaminide -- Nagluc
esculetin plus ions-	(substantially black)	N-acetyl + other sugar components
		butyrate -- esterase
		caprylate -- esterase
		palmitate -- esterase

Table VI is a concise summary of the specific examples.

TABLE VI

EXAMPLE SUMMARIES					
Example #	Substrate	<i>E. coli</i>	General Coliforms	<i>Aeromonas</i>	<i>Salmonella/Shigella</i>
I	8-hydroxyquinoline- β -D-glucuronide	X			
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	color \Rightarrow	Black	Purple-Blue	Pink	Teal
II	8-hydroxyquinoline- β -D-glucuronide	X			
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	6-chloro-3-indolyl- β -D-mannoside	X	X		X
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
IIIA	color \Rightarrow	Black	Purple-blue	Red-pink	Purple-blue
	8-hydroxyquinoline- β -D-glucuronide	X			
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	6-chloro-3-indolyl- α -D-galactopyranoside	X	X		X

TABLE VI-continued

EXAMPLE SUMMARIES					
Example #	Substrate	<i>E. coli</i>	General Coliforms	<i>Aeromonas</i>	<i>Salmonella/Shigella</i>
IIIB	5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	color \Rightarrow	Black	Purple-blue	Purple-blue	Pink
	8-hydroxyquinoline- β -D-glucuronide	X			
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
IIIC	5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	color \Rightarrow	Black	Purple-blue	Purple-blue	Teal
May eliminate <i>Aeromonas</i> with inhibitors which allows removal of 6-chloro-3-indolyl- β -D-galactopyranoside from Examples IIIA and IIIB					
IV	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	6-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	color \Rightarrow	Purple-blue	Purple-blue	Pink	Purple-blue
V-A	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	6-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	color \Rightarrow	Purple-blue	Purple-blue	Purple-blue	Pink
V-B	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	color \Rightarrow	Purple-blue	Purple-blue	Purple-blue	Teal
V-C May eliminate <i>Aeromonas</i> with inhibitors which allows removal of substrate No. 1 from example V-A and allows removal of substrate No. 3 from example V-B					
VI	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	color \Rightarrow	Purple-blue	Purple-blue	Pink	Teal
VII	8-hydroxyquinoline- β -D-glucuronide	X			
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	6-chloro-3-indolyl-N-acetyl- β -D-galactosaminide		X	X	
	Note: In example 7, <i>Vibrio</i> and <i>Plesiomonas</i> also show as pink along with <i>Aeromonas</i>				
	color \Rightarrow	Black	Purple-blue	Pink (see note)	Teal
VIII	8-hydroxyquinoline- β -D-glucuronide	X			
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X		X
	6-chloro-3-indolyl-N-acetyl- β -D-galactosaminide		X	X	
	Note: In example 8, <i>Vibrio</i> and <i>Plesiomonas</i> also show as pink along with <i>Aeromonas</i>				
IX	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	6-chloro-3-indolyl-N-acetyl- β -D-galactopyranoside	X	X	X	
	color \Rightarrow	Purple-blue	Purple-blue	Pink	Teal
X	6-chloro-3-indolyl-N-acetyl- β -D-galactopyranoside		X	X	
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	Note: For example 10, <i>Vibrio</i> and <i>Plesiomonas</i> also show as pink along with <i>Aeromonas</i>				
XI	8-hydroxyquinoline- β -D-glucuronide	X			
	6-chloro-3-indolyl- β -D-galactopyranoside (or)	X	X	X	
	5-bromo-6-chloro-3-indolyl- β -D-galactopyranoside	X	X	X	
	Note: In example 11, <i>Aeromonas</i> may be eliminated by adding inhibitors				
	color \Rightarrow	Black	Pink or Teal	Pink or Teal	Not detected

TABLE VI-continued

EXAMPLE SUMMARIES					
Example #	Substrate	<i>E. coli</i>	General Coliforms	<i>Aeromonas</i>	<i>Salmonella/Shigella</i>
XII	8-hydroxyquinoline- β -D-galactopyranoside	X	X	X	
	5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	(or) 6-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
	Note: In example 13, <i>Aeromonas</i> may be eliminated by adding an inhibitor				
XIII	Use same substrates as in example No. 1, and add: 4-methyl-umbelliferyl- β -D-xylopyranoside	Black Enterobacter and Klebsiella showing as black colonies will fluoresce, thereby allowing reduction in false positive count of <i>E. coli</i> .	Black	Black	Teal or Pink
XIV	8-hydroxy-quinoline- β -D-glucuronide	X			
XV	6-chloro-3-indolyl-caprylate				X Red-pink
	8-hydroxy-quinoline- β -D-glucuronide	Black			
	5-bromo-4-chloro-3-indolyl-caprylate	X			X
	6-chloro-3-indolyl- α -D-galactopyranoside	X	X		X
XVI	8-hydroxy-quinoline- β -D-glucuronide	Black	Red-pink		Blue-violet
	5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside	X	X	Inhibitor present	
	(or)	X	X		
	5-bromo-6-chloro-3-indolyl- β -D-galactopyranoside	X	X		
	color \Rightarrow	Black	Teal (or) Magenta		

Example I

The microorganisms chosen to be identified, quantified and differentiated are *E. coli* general coliforms, *Aeromonas*, *Shigella* or *Salmonella*.

With reference to Table IV, *E. coli* produces the enzyme Bgluc, and Bgluc is not produced by any of the other microorganisms desired to be detected. With reference to the right hand side of Table V, it can be seen that the test enzyme Bgluc has a corresponding substrate component of β -D-glucuronide. Thus, a chromogenic or nonchromogenic component which produces a distinct color upon cleavage of Bgluc should be chosen from the left hand side of Table V. 8-hydroxyquinoline is chosen for its preferred substantially black color. The first identified substrate is therefore 8-hydroxyquinoline- β -D-glucuronide, the availability of which is described above. A metallic salt such as ferric ammonium citrate is also required and is added to the test medium so that, upon cleavage of the substrate by Bgluc, a substantially black water-insoluble complex forms in the medium. The substantially black precipitate consists of the ferric ions and the aglycone released when the substrate is hydrolyzed by the glucuronidase from *E. coli*.

With further reference to Table IV, Bgal, Bfuc and Bglu are common to *Aeromonas* and general coliforms. However, as indicated in Table IV, Bgal, Bfuc and Bglu are not produced generally by *Salmonella* and *Shigella*. Therefore, a substrate component corresponding to one of Bgal, Bfuc and Bglu can be selected from the right hand side of Table V. Bgal and the associated substrate component β -D-galactopyranoside are chosen. The 6-chloro-3-indolyl-chromogenic component produces a red-pink color upon cleavage from its substrate in the presence of Bgal and is selected as the chromogenic component. The second substrate is therefore 6-chloro-3-indolyl- β -D-galactopyranoside.

Again referring to Table IV, Bman, Aara and Agal are common to *Salmonella*, *Shigella* and general coliforms.

However, as indicated in Table IV, Bman, Aara and Agal are not produced by *Aeromonas*. Thus, one of Bman, Aara and Agal can be chosen and its associated substrate component identified with reference to Table V. The test enzyme Agal and the respective substrate component α -D-galactopyranoside are chosen. Next, a chromogenic component must be selected from Table V. As shown on the left hand side of Table V, the chromogenic component 5-bromo-4-chloro-3-indolyl produces a teal color upon cleavage from its associated substrate and is therefore selected. The third substrate is therefore 5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside.

General coliforms have a wide enzyme pattern which is responsive to both the 6-chloro-3-indolyl- β -D-galactopyranoside substrate and the 5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside substrate. Therefore, general coliforms will show as a fourth distinct color which is a combination of the colors produced by the two aforementioned substrates, respectively. In this case the fourth color will be violet-blue, which is a combination of red-pink and teal.

Finally, as seen in Table IV, *E. coli* also exhibits a wide enzyme pattern and responsive to all three of the substrates chosen in this example, viz., 8-hydroxyquinoline- β -D-glucuronide, 6-chloro-3-indolyl- β -D-galactopyranoside, and 5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside. Nonetheless, *E. coli* colonies present in the test medium will show as a substantially black color because, as discussed above, the chromogenic substrates do not interfere with the substantially black color. Advantageously, this substantially black color provides a superior means for distinguishing the *E. coli*, as well as allows four separate microorganisms to be detected, quantified, differentiated and identified in a single test medium. See Table VI.

Example II

The selected microorganisms to be detected, quantified, differentiated and identified are *E. coli* as a first color;

general coliforms, *Salmonella* and *Shigella* as a second color; and *Aeromonas* as a third color.

With reference to Table IV, *E. coli* produces the enzyme Bgluc, and Bgluc is not produced by any of the other microorganisms desired to be detected. With reference to the right hand side of Table V, it can be seen that the test enzyme Bgluc has a corresponding substrate component of β -D-glucuronide. Thus, a chromogenic or nonchromogenic component which produces a distinct color upon cleavage of Bgluc should be chosen from the left hand side of Table V. 8-hydroxyquinoline is chosen for its preferred substantially black color. The first identified substrate is therefore 8-hydroxyquinoline- β -D-glucuronide, the availability of which is described above. A metallic salt such as ferric ammonium citrate is also required and is added to the test medium so that, upon cleavage of the substrate by Bgluc, a substantially black water-insoluble complex forms in the medium. The substantially black precipitate consists of the ferric ions and the aglycone released when the substrate is hydrolyzed by the glucuronidase from *E. coli*.

With further reference to Table IV, Bgal, Bfuc and Bglu are common to *Aeromonas* and general coliforms. However, as indicated in Table IV, Bgal, Bfuc and Bglu are not produced by *Salmonella* or *Shigella*. Using Table V in the fashion described above, 6-Chloro-3-indolyl- β -D-galactopyranoside is selected as the second substrate, which will produce a red-pink color upon cleavage as indicated by the chromogenic component list of Table V.

As seen in Table IV, the enzyme Bman is common to *Salmonella/Shigella* but not *Aeromonas*. From table V, the substrate component associated with Bman is β -D-mannopyranoside. In this example, it is desired to also produce the second distinct color (red-pink) with *Salmonella/Shigella* so that, ultimately, *Salmonella/Shigella* colonies present in the test medium will show as the same color as general coliforms present in the test medium. Thus, the chromogenic component is 6-Chloro-3-indolyl- and the third substrate is therefore 6-Chloro-3-indolyl- β -D-mannopyranoside.

In this example, again using Table V, a fourth substrate is identified that will be cleaved by one of the enzymes Bman, Aara, Agal common to *Salmonella/Shigella* to produce a third distinct color. Using table V in the fashion described above, the fourth substrate selected is 5-Bromo-4-chloro-3-indolyl- α -D-galactopyranoside, which produces a teal-green color in the presence of Agal common to *Salmonella/Shigella*.

The resulting colors of colonies present in the test medium can be predicted as follows. *E. coli* exhibits a wide enzyme pattern that is positive for all four of the substrates chosen in this example, including the 8-hydroxy-glucuronide substrate which produces a substantially black color upon cleavage in the presence of the ions of the ferric salt. *E. coli* colonies show as substantially black. *Aeromonas* has an enzyme pattern which reacts with only the 6-Chloro-3-indolyl- β -D-galactopyranoside substrate chosen in this example and therefore colonies of *Aeromonas* show as red-pink. *Salmonella/Shigella* has an enzyme pattern which cleaves both the third and fourth substrates selected in this example and therefore colonies of *Salmonella/Shigella* show as purple-blue (a combination of teal and red-pink). Finally,

general coliforms are positive for each of the second, third and fourth substrates selected and colonies thereof show as purple-blue, indistinguishable from the *Salmonella/Shigella* colonies. As discussed above, different strains of all species of the various genera will not all produce the same amounts of the various enzymes, so there may be slight variations in shades of purple-blue, for example.

Example IIIA

The selected microorganisms to be quantified and differentiated in this example are *E. coli* as a first color, general coliforms and *Aeromonas* as a second color, and *Salmonella/Shigella* as a third color.

With reference to Table IV, *E. coli* produces the enzyme Bgluc, and Bgluc is not produced by any of the other microorganisms desired to be detected. With reference to the right hand side of Table V, it can be seen that the test enzyme Bgluc has a corresponding substrate component of β -D-glucuronide. Thus, a chromogenic or nonchromogenic component which produces a distinct color upon cleavage of Bgluc should be chosen from the left hand side of Table V. 8-hydroxyquinoline is chosen for its preferred substantially black color. The first identified substrate is therefore 8-hydroxyquinoline- β -D-glucuronide, the availability of which is described above. A metallic salt such as ferric ammonium citrate is also required and is added to the test medium so that, upon cleavage of the substrate by Bgluc, a substantially black water-insoluble complex forms in the medium. The substantially black precipitate consists of the ferric ions and the aglycone released when the substrate is hydrolyzed by the glucuronidase from *E. coli*.

Using tables IV and V in a fashion similar to that described above with reference to Examples I and II, 6-Chloro-3-indolyl- β -D-galactopyranoside is selected as a second substrate to combine with one of the enzymes Bgal, Bfuc and Bglu common to coliforms and *Aeromonas*, but negative for *Salmonella/Shigella* to produce a second distinct color, in this case substantially red-pink.

Similarly, 6-Chloro-3-indolyl- α -D-galactopyranoside is selected as a third substrate to combine with Agal, which is common to coliforms and *Salmonella/Shigella*, but negative for *Aeromonas*. Upon reaction with the enzyme, this substrate will also produce the same distinct second color, namely red-pink.

5-Bromo-4-chloro-3-indolyl- β -D-galactopyranoside is selected as a fourth substrate to combine with the enzyme Bgal, which is common to coliforms and *Aeromonas*, but negative for *Salmonella/Shigella*. This fourth substrate produces a teal-green color upon reaction with Bgal.

The resulting colors of colonies present in the test medium can be predicted as follows. *E. coli* exhibits a wide enzyme pattern and is positive for all four of the substrates chosen in this example. Therefore, *E. coli* colonies will show as substantially black. General coliform colonies have an enzyme pattern which is positive for the second, third and fourth substrates, so that general coliform colonies show as purple-blue. *Aeromonas* colonies have an enzyme pattern which is positive for the second and fourth substrates chosen so that *Aeromonas* colonies also show as purple-blue. Finally, the enzymes common to *Salmonella/Shigella* are

only positive for the third of the four substrates, so that Salmonella/Shigella colonies show as red-pink.

Example IIIB

As a variation, the test medium of Example IIIA can be prepared such that colonies of Salmonella/Shigella will show as teal instead of pink-red, all of the other colony colors being the same as Example IIIA. With reference to Table VI, such can be accomplished by replacing the 6-chloro-3-indolyl- α -D-galactopyranoside of Example IIIA with 5-bromo-4-chloro-3-indolyl- α -D-galactopyranoside.

Example IIIC

A second, independent method for producing the same three colors as Example IIIA for the same four components can be achieved by adding nalidixic acid or other antibiotics or inhibitors of Aeromonas to the components listed in Example I. In so doing, the cefsulodin or nalidixic acid or other substance acts as an inhibitor for Aeromonas so Aeromonas colonies do not grow. If Aeromonas is eliminated, then the purple-blue colonies are all true coliforms. If not eliminated, any Aeromonas will be counted as part of the coliforms which some persons may prefer since Aeromonas is an important indicator organism.

Example IV

In this example, the selected microorganisms to be detected, quantified and differentiated are *E. coli* coliforms and Salmonella/Shigella as a first distinct color and Aeromonas as a second distinct color. One test medium which achieves this result is the test medium described in Example II, except that the first substrate and metallic salt are omitted. Thus, because the enzyme pattern of *E. coli* reacts with the same substrates as the enzyme pattern for general coliforms, *E. coli* and general coliforms will be the same color in this test medium. Specifically, *E. coli*, coliforms and Salmonella/Shigella colonies will show as a purple-blue color, whereas Aeromonas colonies will show as a substantially red-pink color.

Example V

The selected microorganisms to be detected, quantified and differentiated are *E. coli*, general coliforms and Aeromonas as a first distinct color, and Salmonella/Shigella as a second distinct color. One test medium which achieves this result is the test medium of Example 3 with the first substrate and metallic salt being omitted. In this test medium, *E. coli*, general coliforms and Aeromonas colonies will show as a generally purple-blue color, whereas Salmonella and Shigella colonies will show as a generally teal-green color or as a red-pink color.

Optionally, the 6-chloro-3-indolyl- α -D-galactoside can be replaced with 5-bromo-4-chloro-3-indolyl- β -D-galactoside so that Salmonella colonies show as teal, rather than pink.

A third way to achieve the same result is with an antibiotic, preferably nalidixic acid, to inhibit the growth of Aeromonas colonies. If Aeromonas is eliminated, then the purple-blue colonies are all true coliforms. If not eliminated, any Aeromonas will be counted as part of the coliforms

which some persons may prefer since Aeromonas is an important indicator organism.

Example VI

The selected microorganisms to be detected, quantified and differentiated are *E. coli* and coliforms as a first distinct color, Aeromonas as a second distinct color and Salmonella/Shigella as a third distinct color. A test medium which achieves this result is the test medium of Example I with the first substrate and metallic salt being omitted. In such a test medium, *E. coli* and general coliform colonies will show as purple-blue, Aeromonas colonies will show as generally red-pink, and Salmonella or Shigella colonies will show as generally teal-green.

Example VII

The selected microorganisms to be detected, quantified and differentiated are *E. coli* as a first distinct color which is substantially black; general coliforms as a second distinct color which is substantially purple-blue; Aeromonas/Vibrio/Plesiomonas as a third distinct color which is substantially red-pink; and Salmonella or Shigella as a fourth distinct color which is substantially teal-green.

With reference to Table IV, *E. coli* produces the enzyme Bgluc, and Bgluc is not produced by any of the other microorganisms desired to be detected. Therefore, a substrate which produces a distinct color upon cleavage of Bgluc should be chosen from Table V. 8-hydroxyquinoline- β -D-glucuronide produces a substantially black color in the presence of Bgluc and would be the preferred choice of substrate, as explained below. A metallic salt such as ferric ammonium citrate is also added to form a substantially black water insoluble complex consisting of the ferric ions and the aglycone released when the substrate is hydrolyzed by the glucuronidase from *E. coli*.

With further reference to Table IV, it can be seen that the enzyme Ngai and Naglu are common to the microorganisms Aeromonas, Plesiomonas, and Vibrios. Therefore, a suitable substrate for testing all of these microorganisms as a single distinct color is 6-chloro-3-indolyl-N-acetyl- β -D-galactosaminide, which produces a substantially red-pink color in the presence of these enzymes.

Again referring to Table IV, Bman, Aara and Agal are common to Salmonella/Shigella and general coliforms. However, as indicated in Table IV, Bman, Aara and Agal are not produced by Aeromonas. Therefore, a substrate can be selected from Table V which reacts with one of Bman, Aara and Agal to produce a third distinct color. As shown in Table V, 5-bromo-4-chloro-3-indolyl- α -D-galactoside produces a teal-green color in the presence of Agal and is therefore selected as a substrate.

In this test medium *E. coli* colonies will show as substantially black, general coliform colonies will show as substantially purple-blue, Aeromonas, Vibrio and Plesiomonas colonies will show as substantially red-pink, and Salmonella or Shigella colonies will show as substantially teal.

Example VIII

The selected microorganisms to be detected, quantified and differentiated are *E. coli* as a first distinct color;

33

coliforms, *Salmonella* or *Shigella* as a second distinct color; and *Aeromonas*, *Vibrio* and *Plesiomonas* as a third distinct color. One test medium for achieving this result is the test medium of Example 2, except that the fourth substrate chosen is 6-Chloro-3-indolyl-N-acetyl- α -D-galactosaminide, to which each of the microorganisms *Plesiomonas*, *Vibrios* and *Aeromonas* are responsive so that each of these colonies shows as a generally red-pink color.

Example IX

The selected microorganisms to be detected, quantified and differentiated in this example are *E. coli* and general coliforms as a first distinct color which is purple-blue; *Aeromonas*, *Plesiomonas*, and *Vibrios* as a second distinct color which is red-pink; and *Salmonella* or *Shigella* as a third distinct color which is teal-green. This result can be achieved with the test medium as described in Example 6 with the addition of 6-Chloro-3-indolyl-N-acetyl- β -D-galactosaminide, to which each of the microorganisms *Plesiomonas*, *Vibrio* and *Aeromonas* is responsive.

Example X

The selected microorganisms to be detected, quantified and differentiated in this example are *E. coli* and general coliforms as a first color; *Aeromonas*, *Vibrio* and *Plesiomonas* as a second distinct color; and *Salmonella* or *Shigella* as a third distinct color. A suitable test medium that achieves this result is the test medium disclosed in Example 7 except that the first substrate for detecting *E. coli* colonies is omitted. In this example, *E. coli* and general coliform colonies show as generally purple-blue, *Aeromonas*, *Vibrio* and *Plesiomonas* show as generally red-pink, and *Salmonella* or *Shigella* show as generally teal-green. The addition of 6-Chloro-3-indolyl- β -D-galactopyranoside is necessary to yield the purple-blue color for *E. coli* colonies.

Example XI

The selected microorganisms to be detected, quantified and differentiated in this example are *E. coli* as a substantially black color and general coliforms as a red-pink color.

With reference to Table IV, *E. coli* produces the enzyme Bgluc, and Bgluc is not produced by any of the other microorganisms desired to be detected. Therefore, a substrate which produces a distinct color upon cleavage of Bgluc should be chosen from Table V. 8-hydroxyquinoline- β -D-glucuronide produces a dark color in the presence of Bgluc and would be the preferred choice of substrate. A metallic salt such as ferric ammonium citrate is also added to form a black water insoluble complex consisting of ferric ions and the aglycone released when the substrate is hydrolyzed by glucuronidase from *E. coli*.

With further reference to Table IV, Bgal, Bfuc and Bglu are common to *Aeromonas* and general coliforms. However, as indicated in Table IV, Bgal, Bfuc and Bglu are not generally produced by *Salmonella* or *Shigella*. Therefore, a substrate can be selected from Table V which reacts with one of Bgal, Bfuc and Bglu to produce a second distinct color. 6-chloro-3-indolyl- β -D-galactopyranoside produces a pink color in the presence of Bgal and is selected as the second substrate.

34

Optionally, the 6-chloro-3-indolyl- β -D-galactopyranoside can be replaced with 5-bromo-6-chloro-3-indolyl- β -D-galactopyranoside so that *Aeromonas* and general coliform colonies show as teal instead of pink.

As noted, the second substrate selected will result in colonies of *Aeromonas* also showing as a generally red-pink color. To avoid growth of *Aeromonas* colonies, an inhibitor, preferably nalidixic acid, is added. Thus, colonies of *E. coli* will show as substantially black, whereas colonies of general coliforms will show as a red-pink color.

Example XII

The selected microorganisms to be detected, quantified and differentiated in this example are *E. coli*, general coliforms and *Aeromonas* spp. as a substantially black color and *Salmonella* or *Shigella* spp. as a second distinct color. The first substrate selected is 8-hydroxyquinoline- β -D-galactoside, which results in colonies of *E. coli*, general coliforms and *Aeromonas* showing as substantially black. The second substrate chosen can be either 5-Bromo-4-chloro-3-indolyl- α -D-galactopyranoside or 6-chloro-3-indolyl- α -D-galactopyranoside. If the former of these two substrates is chosen, colonies of *Salmonella* or *Shigella* will show as a teal color, whereas if the latter of the two aforementioned substrates is chosen, colonies of *Salmonella* or *Shigella* will show as a red-pink color.

Optionally, in this example, *Aeromonas* may be eliminated by adding an inhibitor, preferably nalidixic acid, as discussed in detail above.

Example XIII

The selected microorganisms to be detected, quantified and differentiated in this example are the same as in Example 1, except that this example illustrates a correction for false positives. That is, it is possible that certain unusual *Enterobacter* and *Klebsiella* spp. will show as black colonies along with *E. coli* in the test medium disclosed in Example 1. Thus, the count of *E. coli* could be inaccurately high.

In this example, 4-methyl-umbelliferyl- β -D-xylopyranoside is added to the test medium described in Example 1. In so doing, *Enterobacter* and *Klebsiella* spp. showing as black colonies will also fluoresce, thereby allowing reduction in the false positive count of *E. coli*. This example illustrates the flexibility of embodiments incorporating the present invention. The fluorescent component does not interfere with the substantially black color so that the black colonies are easily distinguished with the naked eye. Yet, under ultraviolet light, false positives can be detected and substantially reduced by examining the black colonies for fluorescence.

Example XIV

The selected microorganisms to be detected, quantified, differentiated and identified are *E. coli* as a substantially black color and *Salmonella* spp. as pink-red. General coliforms are colorless in this example.

With reference to Example 1, *E. coli* is responsive to 8-hydroxy-quinoline- β -D-glucuronide. General coliforms, *Salmonella*, *Shigella* and *Aeromonas* are not responsive to 8-hydroxy-quinoline- β -D-glucuronide. Thus, the first substrate chosen is 8-hydroxy-quinoline- β -D-glucuronide.

With reference to table IV, esterase enzyme is positive for *Salmonella* spp., but not any of the other preferred microorganisms to be detected. With reference to table V, the substrate 6-chloro-3-indolyl-caprylate can be identified, and will produce a pink-red color upon cleavage, and is therefore chosen as the second substrate.

In this test medium, colonies of *E. coli* will show as substantially black and colonies of *Salmonella* or *Shigella* will show as pink-red.

Example XV

The selected microorganisms to be detected, quantified, differentiated and identified are *E. coli* as a substantially black color, *Salmonella* spp. as dark blue-purple, and general coliforms as red-pink.

With reference to Example I, *E. coli* is responsive to 8-hydroxy-quinoline- β -D-glucuronide. General coliforms, *Salmonella*, *Shigella* and *Aeromonas* are not responsive to 8-hydroxy-quinoline- β -D-glucuronide. Thus, the first substrate chosen is 8-hydroxy-quinoline- β -D-glucuronide.

With reference to tables IV and V, 5-bromo-4-chloro-3-indolyl-caprylate can be identified as the second substrate to which *Salmonella* or *Shigella* will be responsive. With further reference to Table V, 5-bromo-4-chloro-3-indolyl-caprylate forms a teal color upon cleavage.

6-chloro-3-indolyl- α -D-galactopyranoside, which produces a pink-red color upon cleavage, is chosen as the third substrate to which *E. coli* general coliforms and *Salmonella* or *Shigella* are responsive.

In this example, *E. coli* colonies show as substantially black, general coliform colonies show as red-pink, and *Salmonella*/*Shigella* show as blue-violet (=red-pink+teal).

Example XVI

The selected microorganisms to be detected, quantified and differentiated in this example are *E. coli* as a substantially black color and general coliforms as a second distinct color.

With reference to Table IV, *E. coli* produces the enzyme Bgluc, and Bgluc is not produced by any of the other microorganisms desired to be detected. Therefore, a substrate which produces a distinct color upon cleavage of Bgluc should be chosen from Table V. 8-hydroxyquinoline- β -D-glucuronide produces a dark color in the presence of Bgluc and would be the preferred choice of substrate. A metallic salt such as ferric ammonium citrate is also added to form a black water insoluble complex consisting of ferric ions and the aglycone released when the substrate is hydrolyzed by glucuronidase from *E. coli*.

With further reference to Table IV, Bgal, Bfuc and Bglu are common to *Aeromonas* and general coliforms. However, as indicated in Table IV, Bgal, Bfuc and Bglu are not generally produced by *Salmonella*/*Shigella*. Therefore, a substrate can be selected from Table V which reacts with one of Bgal, Bfuc and Bglu to produce a second distinct color. 5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside can be chosen as the second substrate, in which event colonies of *E. coli* will appear as substantially black and general coliform colonies will appear as teal. Optionally, 5-bromo-6-chloro-3-indolyl-galactopyranoside can be chosen as the second

substrate, in which event colonies of *E. coli* will appear as substantially black and general coliform colonies will appear magenta. To avoid growth of *Aeromonas* colonies, an inhibitor, preferably nalidixic acid, is added. Thus, colonies of *E. coli* will show as substantially black, whereas colonies of general coliforms will show as a magenta color.

Although several broad examples which incorporate the present invention have been described above, it is to be understood that the present invention is not to be limited by the examples disclosed herein. Indeed, the disclosure and examples above teach one of ordinary skill a virtually limitless number of test media which would be within the scope of the claims appended hereto.

Further, while this invention has been described as having a preferred design, the present invention can be further modified within the spirit and scope of this disclosure. This application is therefore intended to cover any variations, uses, or adaptations of the invention using its general principles. Further, this application is intended to cover such departures from the present disclosure as come within known or customary practice in the art to which this invention pertains and which fall within the limits of the appended claims.

What is claimed is:

1. A test medium for detecting, quantifying or differentiating general coliforms, *E. coli*, *Aeromonas* and *Salmonella*, said test medium comprising: a nutrient base medium including ions of a salt; a first substrate which forms a first component of a first color in the presence of *E. coli*; a second substrate which forms a second component of a second color in the presence of *Salmonella*; and a third substrate which forms a third component of a third color in the presence of *Aeromonas*; said second and third substrate forming said second and third components, respectively to make a fourth color, which is a combination of said second and third colors, in the presence of general coliforms, all of said colors being distinguishable from one another; said first substrate being a β -glucuronide nonchromogenic substrate; said second substrate being an α -D-galactoside chromogenic substrate and said third substrate being a β -D-galactoside chromogenic substrate, and said first color being substantially black.

2. The test medium of claim 1, wherein said first substrate is selected from the group consisting of 8-hydroxyquinoline- β -D-glucuronide, an esculetin glucuronide, and cyclohexenoesculetin- β -D-glucuronide.

3. The test medium of claim 1, wherein said first substrate is 8-hydroxyquinoline-3-D-glucuronide and forms a substantially nondiffusible compound in the presence of ions of said salt and *E. coli*.

4. The test medium of claim 1, wherein said salt comprises a metallic salt and said first component is water insoluble as formed by reaction with said ions.

5. The test medium of claim 1, wherein said first substrate consists essentially of 8-hydroxyquinoline- β -D-glucuronide, said second substrate consists essentially of 5-bromo-4-chloro-3-indole- β -D-galactoside, and said third substrate consists essentially of 6-chloro-3-indole- β -D-galactoside.

6. A method for detecting, quantifying, or differentiating colonies of *Aeromonas* from selected other biological enti-

37

ties in a test sample, said method comprising the following steps: providing a base medium including ions of salt, a β -D-galactoside substrate that forms a first component of a first color in the presence of a first enzyme, an α -D-galactoside substrate that forms a second component of a second color distinguishable from said first color in the presence of a second enzyme, and a β -glucuronide nonchromogenic substrate that forms a third substantially black component in the presence of a third enzyme; inoculating the test medium with a test sample; incubating the test medium; and examining the test medium whereby aggregations of colonies of *Aeromonas* are indicated by said first

38

color, and aggregations of colonies of *Salmonella* are indicated by said second color, and whereby colonies of general coliforms are indicated by a third color, said third color being a combination of said first and second colors.

7. The method as set forth in claim 6, further comprising the step of examining the test medium for *E. coli* as indicated by the presence of substantially black aggregates.

8. The method as set forth in claim 7, wherein said β -glucuronide substrate is 8-hydroxyquinoline.

* * * * *

UNITED STATES PATENT AND TRADEMARK OFFICE
CERTIFICATE OF CORRECTION

PATENT NO. : 6,787,332 B2
DATED : September 7, 2004
INVENTOR(S) : Geoffrey N. Roth et al.

Page 1 of 1


It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

Column 36.

Lines 59-65, delete "β" and insert -- α -- (all occurrences)

Signed and Sealed this

Twenty-fourth Day of May, 2005

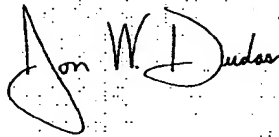A handwritten signature in black ink, appearing to read "Jon W. Dudas", is written over a rectangular area of the document.

JON W. DUDAS
Director of the United States Patent and Trademark Office

**TRANSMITTAL LETTER
(General - Patent Issued)**

Docket No.
MLB0001.01

Patentee(s): Geoffrey N. Roth et al.

U.S. Patent No.
6,787,332

Issue Date
September 7, 2004

Title: Test Media and Quantitative or Qualitative Method for Identification and Differentiation of E. Coli
Coliforms, Salmonella, and Aeromonas in a Test Sample

COMMISSIONER FOR PATENTS:

Transmitted herewith is:

Certificate of Correction

- ☒ No additional fee is required.
- ☐ A check in the amount of _____ is attached.
- ☐ The Director is hereby authorized to charge and credit Deposit Account as described below.
- ☐ Charge the amount of _____
- ☐ Credit any overpayment.
- ☐ Charge any additional fee required.
- ☐ Payment by credit card. Form PTO-2038 is attached.

WARNING: Information on this form may become public. Credit card information should not be included on this form. Provide credit card information and authorization on PTO-2038.


Signature


Daniel Tychonievich, Reg. No. 41,358
Baker & Daniels
205 W. Jefferson Blvd., Suite 250
South Bend, IN 46601
Tel: (574)234-4149
Fax: (574)239-1900

Dated: 3/2/05

I hereby certify that this correspondence is being deposited with the United States Postal Service with sufficient postage as first class mail in an envelope addressed to "Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450" [37 CFR 1.8(a)] on

3/2/05

(Date)

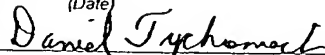

Signature of Person Mailing Correspondence

Daniel Tychonievich

Typed or Printed Name of Person Mailing Correspondence

Customer No.: 27817

CC:

UNITED STATES PATENT AND TRADEMARK OFFICE CERTIFICATE OF CORRECTION

PATENT NO : 6,787,332
 DATED : September 7, 2004
 INVENTOR(S) : Geoffrey N. Roth et al.

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

Page 36, Claim 5, Line 62, delete " β " and insert – α –

Claim should read:

5. The test medium of claim 1, wherein said first substrate consists essentially of 8-hydroxyquinoline- β -D-glucuronide, said second substrate consists essentially of 5-bromo-4-chloro-3-indole- α -D-galactoside, and said third substrate consists essentially of 6-chloro-3-indole- β -D-galactoside.

MAILING ADDRESS OF SENDER:

Eric J. Groen
 Baker & Daniels
 205 West Jefferson Blvd., Suite 250
 South Bend, IN 46601

PATENT NO. 6,787,332

No. of additional copies



This collection of information is required by 37 CFR 1.322, 1.323, and 1.324. The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 1.0 hour to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, P.O. Box 1450, Alexandria, VA 22313-1450. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2.

The U. S. Patent and Trademark Office date-stamp hereon is acknowledgment that on the date stamped, the U. S. Patent and Trademark Office received:

TYPE OF PAPER : Certificate of Correction

APPLICANT(S) : Geoffrey N. Roth *et al.*
PATENT NO. : 6,787,332
ISSUED : September 7, 2004

TITLE : Test Media and Quantitative or Qualitative Method for
Identification and Differentiation of E. Coli, General Coliforms,
Salmonella, and Aeromonas in a Test Sample

ENCLOSURES : Transmittal Form;
Certificate of Correction

DATE MAILED : 03/02/2005
OUR REF. : MLB0001.01 (65514.1)
DT/mr



BDD801 3227327v1

UNITED STATES PATENT AND TRADEMARK OFFICE
CERTIFICATE OF CORRECTION

PATENT NO. : 6,787,332 B2
DATED : September 7, 2004
INVENTOR(S) : Geoffrey N. Roth et al.

Page 1 of 1

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

Column 36,

Lines 59-65, delete "β" and insert -- α -- (all occurrences)


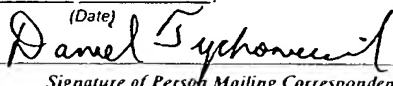


Signed and Sealed this

Twenty-fourth Day of May, 2005

A handwritten signature in cursive script, reading "Jon W. Dudas".

JON W. DUDAS
Director of the United States Patent and Trademark Office

TRANSMITTAL LETTER (General - Patent Issued)		Docket No. MLB0001.01
Patentee(s): Geoffrey N. Roth et al.		
U.S. Patent No. 6,787,332	Issue Date September 7, 2004	
Title: Test Media and Quantitative or Qualitative Method for Identification and Differentiation of E. Coli Coliforms, Salmonella, and Aeromonas in a Test Sample		
<p style="text-align: center;"><u>COMMISSIONER FOR PATENTS:</u></p> <p>Transmitted herewith is:</p> <p style="padding-left: 40px;">Certificate of Correction</p> <p><input checked="" type="checkbox"/> No additional fee is required.</p> <p><input type="checkbox"/> A check in the amount of _____ is attached.</p> <p><input checked="" type="checkbox"/> The Director is hereby authorized to charge and credit Deposit Account 02-0387 as described below.</p> <p style="padding-left: 40px;"> <input type="checkbox"/> Charge the amount of _____ <input checked="" type="checkbox"/> Credit any overpayment. <input checked="" type="checkbox"/> Charge any additional fee required. </p> <p><input type="checkbox"/> Payment by credit card. Form PTO-2038 is attached.</p> <p>WARNING: Information on this form may become public. Credit card information should not be included on this form. Provide credit card information and authorization on PTO-2038.</p> <div style="display: flex; justify-content: space-between;"> <div style="width: 45%;">  <p style="text-align: center; margin-top: -10px;"><i>Signature</i></p> <p>Daniel Tychonievich, Reg. No. 41,358 Baker & Daniels LLP 205 W. Jefferson Blvd., Suite 250 South Bend, IN 46601 Tel: (574)234-4149 Fax: (574)239-1900</p> </div> <div style="width: 45%; text-align: right;"> Dated: 06/10/2005 </div> </div> <div style="display: flex; justify-content: space-between;"> <div style="width: 45%;"> Customer No.: 27817 cc: </div> <div style="width: 45%; border: 1px solid black; padding: 5px;"> <p>I hereby certify that this correspondence is being deposited with the United States Postal Service with sufficient postage as first class mail in an envelope addressed to "Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450" [37 CFR 1.8(a)] on</p> <p style="text-align: center;">06/10/2005 (Date)</p> <p style="text-align: center;"> <i>Signature of Person Mailing Correspondence</i></p> <p style="text-align: center;">Daniel Tychonievich <i>Typed or Printed Name of Person Mailing Correspondence</i></p> </div> </div>		

UNITED STATES PATENT AND TRADEMARK OFFICE CERTIFICATE OF CORRECTION

PATENT NO. : 6,787,332 B2
DATED : September 7, 2004
INVENTOR(S) : Geoffrey N. Roth et al.

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below.

Page 36, Claim 3, Line 53, delete "3" and insert - β -

Page 36, Claim 5, Line 60, delete "a" and insert - β -
Claim 5, Line 62, delete "b" and insert - α -
Claim 5, Line 63, delete "a" and insert - β -

Claim 5 should read:

5. The test medium of claim 1, wherein said first substrate consists essentially of 8-hydroxyquinoline- β -D-glucuronide, said second substrate consists essentially of 5-bromo-4-chloro-3-indole- α -D-galactoside, and said third substrate consists essentially of 6-chloro-3-indole- β -D-galactoside.

MAILING ADDRESS OF SENDER:

Daniel Tychonievich, Baker & Daniels LLP
205 W. Jefferson Blvd., Suite 250, South Bend, IN 46601
Customer Number: 27187

PATENT NO. 6,787,332 B2

No. of additional copies

This collection of information is required by 37 CFR 1.322, 1.323, and 1.324. The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 10 hour to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, P.O. Box 1450, Alexandria, VA 22313-1450. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Attention Certificate of Corrections Branch, Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2.

SBIMANI 193177v1



The U. S. Patent and Trademark Office date-stamp hereon is acknowledgment that on the date stamped, the U. S. Patent and Trademark Office received:

TYPE OF PAPER	:	Certificate of Correction
APPLICANT(S)	:	Geoffrey N. Roth
PATENT NO.	:	6,787,332 B
DATED	:	09/07/2004
TITLE	:	Test Media and Quantitative or Qualitative Method for Identification and Differentiation of Biological Materials in a Test Sample
ENCLOSURE	:	Transmittal Letter; Certificate of Correction
DATE MAILED	:	06/16/2005
OUR REF.	:	MLB0001.01 (65514.1)
DT/mr	:	



BDD801 3227327v1



UNITED STATES DEPARTMENT OF COMMERCE
Patent and Trademark Office
ASSISTANT SECRETARY OF COMMERCE AND
COMMISSIONER OF PATENTS AND TRADEMARKS
Washington, DC 20231

7/26/05
Patent No. : 6787335B2
Inventor(s) : Geoffrey N. Roth
Issued : 9/7/04
Title : TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD FOR
IDENTIFICATION AND DIFFERENTIATION OF BIOLOGICAL MATERIALS IN A
TEST SAMPLE
Atty.doc./File No.

Request for Certificates of Correction

Consideration has been given to your request for the issuance of a Certificate of Correction, for the above – identified patent under the provisions of CFR 1.322.

Inspection of the application for the patent reveals that claim 3 & 5 is printed in accordance with the record. Therefore being no fault on the Patent and Trademark Office, It has no authority to issue a certificate of correction under the provision of 1.322.

In view of the forgoing, your request in this matter, is hereby denied.

Future written correspondence concerning this matter should be filed and directed to Decisions & Certificates of Correction Branch.

Henry Randall
Cecelia Newman
Decisions & Certificates
of Correction Branch
(703) 308-9390 Ext. 108

Daniel Tychonievich, Baker & Daniels LLP
205 w. Jefferson Blvd., Suite 250
South Bend, IN 46601

RECEIVED
JUL 29 2005
BAKER & DANIELS

HR/CBN

BAKER & DANIELS LLP

Operating: South Bend
First Bank Building
205 W. Jefferson, Suite 250
South Bend, Indiana 46601
Phone 574-234-4149

1st SOURCE BANK
SOUTH BEND, INDIANA 46634

CHECK NO. 104310
CHECK DATE
09-23-05

ONE HUNDRED AND 00/100 Dollars

CHECK AMOUNT
\$100.00*****

PAY
TO THE
ORDER OF

Commissioner for Patents

Robert H. Davis
VOID AFTER 180 DAYS

⑈104310⑈ ⑆07121228⑆ 119 721⑈9⑈



The U. S. Patent and Trademark Office date-stamp hereon is acknowledgment that on the date stamped, the U. S. Patent and Trademark Office received:

TYPE OF PAPER : Request for Correction of Patent Certificate

IN RE PATENT OF : Geoffrey N. Roth et al.

PATENT NO. : 6,787,332
ISSUED : 09/07/2004

TITLE : TEST MEDIA AND QUANTITATIVE OR QUALITATIVE METHOD
FOR IDENTIFICATION AND DIFFERENTIATION
OF E. COLI, GENERAL COLIFORMS, SALMONELLA, AND
AEROMONAS IN A TEST SAMPLE

ENCLOSURES

: Request for Correction of Patent Certificate;
Check No. 104310 in the amount of \$100.00;
Exhibits A - G

DATE MAILED : 09/23/2005
OUR REF. : MLB0001.01 (65514.1)
DT/MLWH/mr
BD0801 3227327v1